

Table of Contents

iMotions Lab > Analysis > Analysis (General)

Aggregate Replay	2
------------------------	---

Aggregate Replay


This Help Center article explains how to use the Aggregate replay editor

The Aggregate replay editor is similar to the Replay editor, however, instead of replaying the data from *individual* respondents, the Aggregate replay editor replays the *aggregated* data from all respondents included in the analysis. It is the place where you review the sensor data of a group of respondents, typically with the purpose of spotting group-level patterns in the data. For instance, it could be that you would be interested in knowing which scenes of a film evoke the largest emotional responses across all respondents included in the analysis. It could also be that you want to understand whether there are salient differences in how professional musicians (experts) and amateur musicians (novices) read sheet music? The Aggregate replay editor is further also the place where you generate heatmaps (see also: [Heatmaps](#)).

Working within the Aggregate replay editor is almost identical as working within the Replay editor. Hence, for a comprehensive overview of the functionalities, we mostly refer to the Help Center article about the Replay editor (see also: [Replay](#)). However, there are also few important differences to be aware of.

How to access the Aggregate replay editor

Unlike the Replay editor, the Aggregate replay editor is not available at the 'Study level', but only at the 'Analysis level'. This is the part of the software you enter by creating an analysis, or by clicking the name of an existing analysis in the Library.

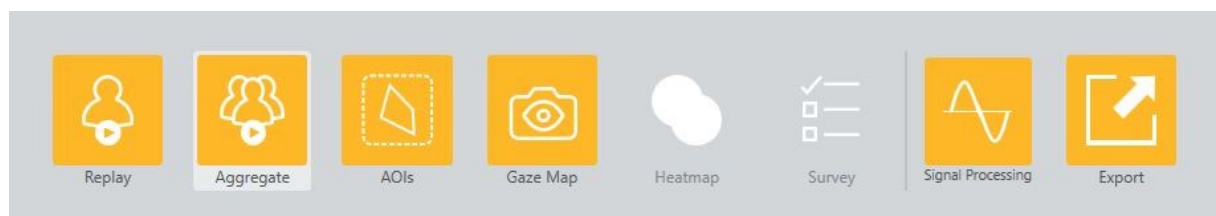


Analysis menu ('Analysis level')

This is because the Replay editor replays the recording of *individual* respondents (consistent with how the data was originally recorded), while the Aggregated replay editor replays *the aggregated data of all respondents included in the analysis*. And thus, in order to compute those aggregations, the analysis must first have been created.

Aggregate replay

Once an analysis is created, the Aggregate replay can be accessed in the following way.

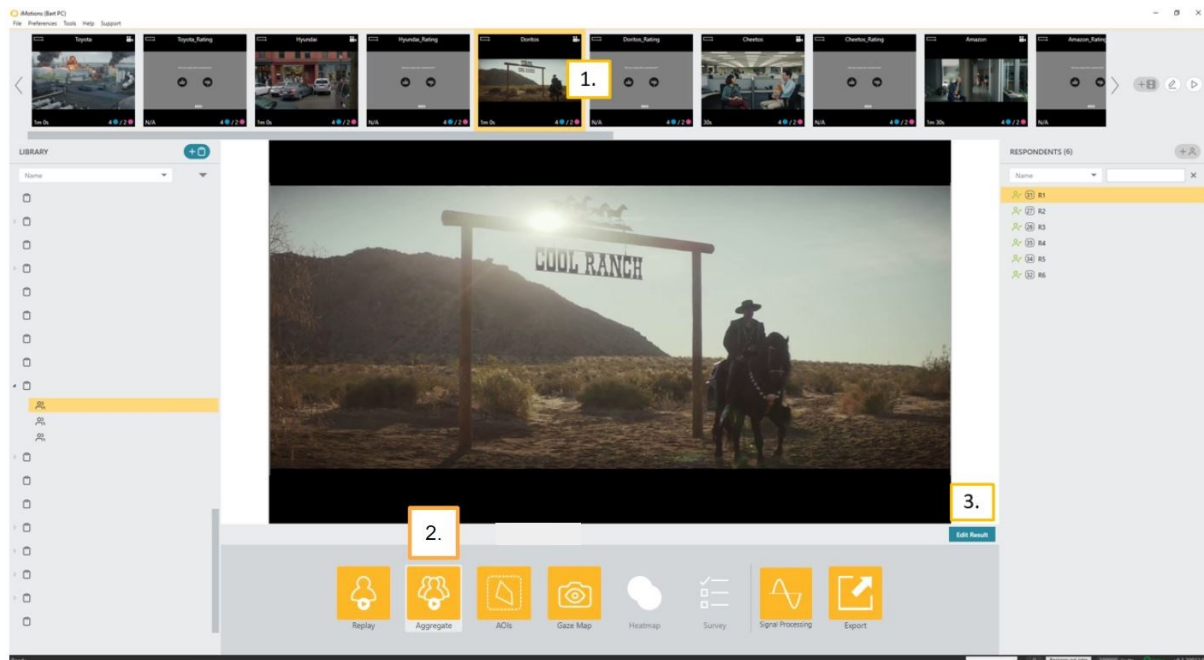

How to access the Aggregate replay editor


Access the Aggregate replay editor (yellow path):

1. Click on a stimulus in the Stimulus bar
2. Click on the option 'Aggregate' in the Analysis bar
3. Click 'Edit result'

Note that the Aggregate replay editor is not available for website recordings, face recordings, scene recordings and screen recordings. The editor is, however, available for 'image scenes', 'video scenes' and 'website scenes'. In the case of a website recording, scene recording or screen recording (but not face recording), it would be possible to first use the scene editor and then access the Aggregate replay editor for the generated 'image scenes', 'video scenes' and 'website scenes' (see also: [Scene editor](#)).


Aggregate replay menu

The menu in the top-left corner of the Aggregate replay editor is nearly identical to the menu in the Replay editor, except that (1) the 'Respondents' menu is not longer present, and that only 'Stimulus annotations' are available.

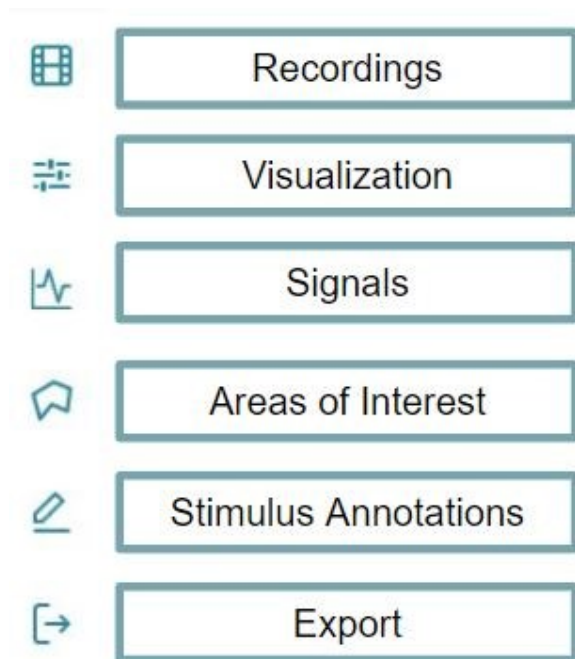

Interface options

The interface options menu bar allows you to customize the Aggregate replay editor, by (un)tickling panels, and by choosing a desired layout. Note that, at the Aggregate replay editor, there are no 'Respondent' and 'Environment' panels; the only available panels are the 'Stimulus', 'Signals', and 'Annotations' panels.

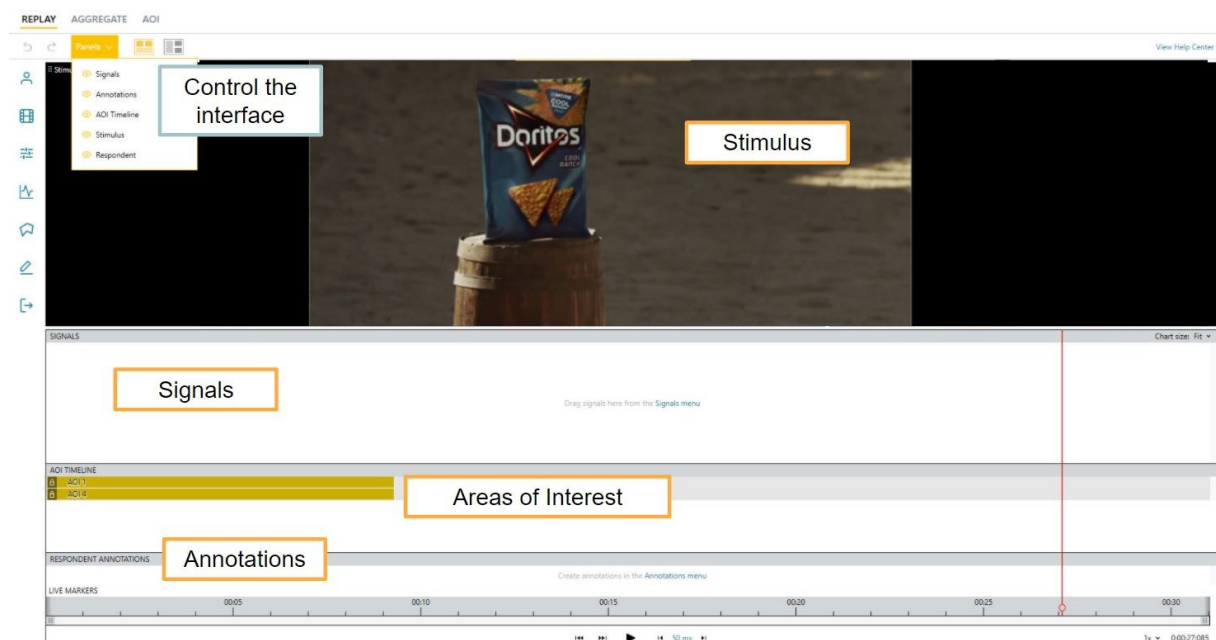


How to control the interface

Control the interface:

1. Navigate to the interface options menu bar
2. (Un)tick the panels of interest
3. Choose a layout

Recordings

The 'Recordings' menu is similar to the 'Recordings' menu in the Replay editor, and allows you to navigate from stimulus to stimulus.

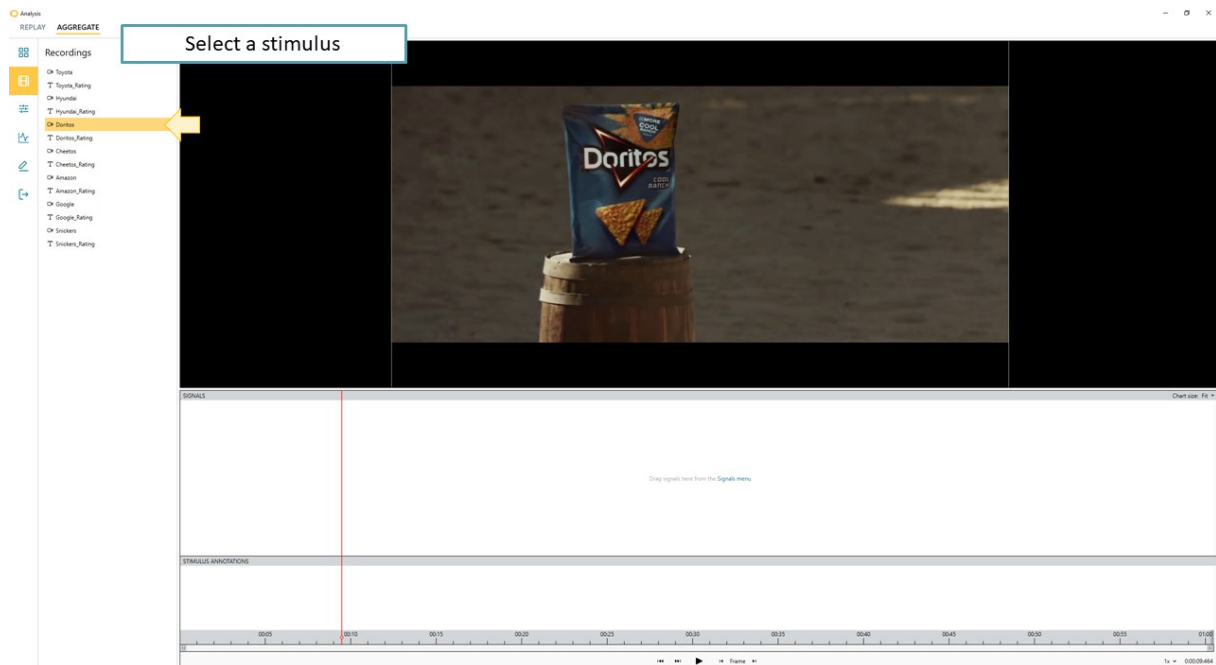


How to navigate from stimulus to stimulus

Select a stimulus:

1. Navigate to the Recordings menu
2. Click on a stimulus in the list

Visualization

In the Replay editor, the 'Visualization' menu is where you generate scan-paths (see also: [Scan-paths](#)); in the Aggregate replay editor, the 'Visualization' menu is where you generate heatmaps. These heatmaps can be static (i.e., summarizing the eye tracking data as an image) or dynamic (i.e., replaying the eye tracking data as a video), and overlay the stimulus in the Stimulus panel (see also: [Heatmaps](#)).

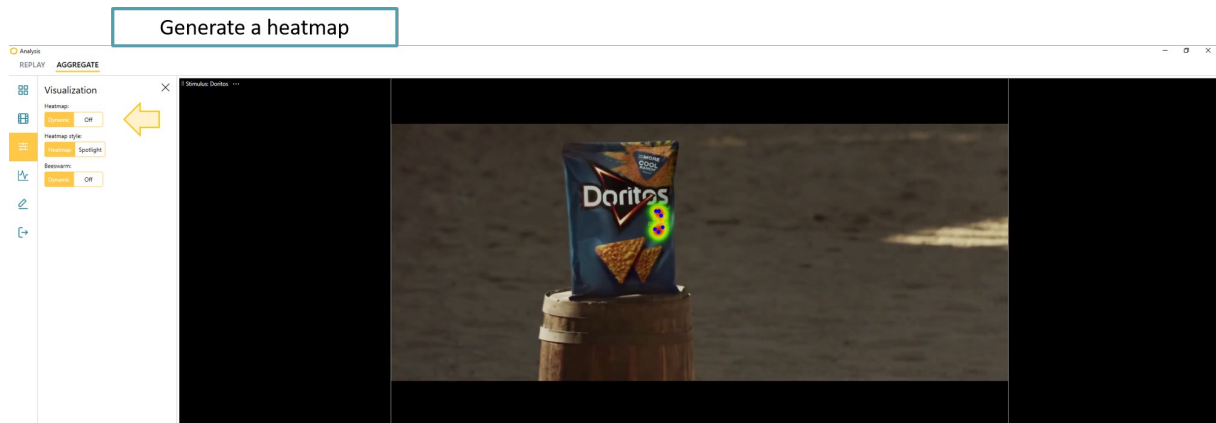

How to generate a heatmap


Generate a heatmap:

1. Navigate to the Visualization menu
2. Consider the options (see also: [Heatmaps](#)).

It is not possible to visualize mouse clicks from within the Aggregate replay editor.

Signals


The 'Signals' menu is again very similar to the Signals menu in the Replay editor, however, in the Aggregated replay editor only the *aggregated* signals are available. These aggregated signals are computed by R-Notebooks, and thus, they will only be available if you have already run the relevant R-Notebooks (see also: [R-Notebooks](#)).

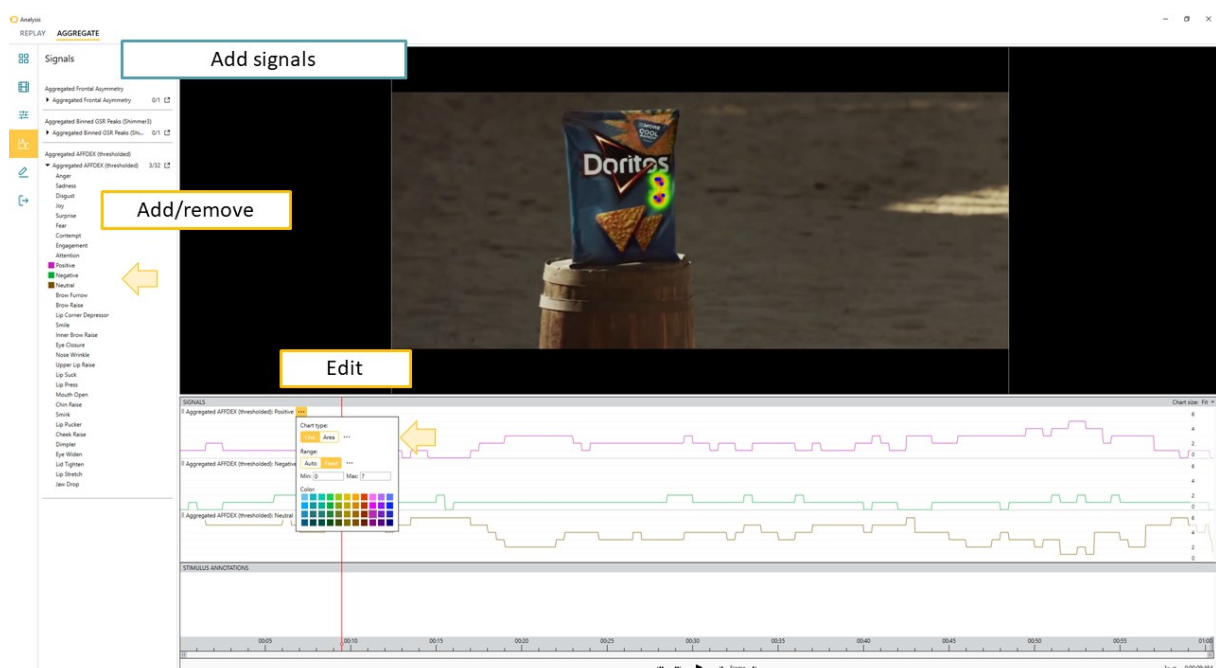

For a comprehensive overview on (1) how to add/remove signals to the Signals panel, and (2) how to edit the graphs, we refer to the Help Center article about the Replay editor (see also: [Replay](#)).


Areas of Interest

The 'Areas of Interest' section allows you to only view what AOIs have been made, and they can only be edited in the [AOI tab](#). Additionally, they can be selected or deselected in order to be shown on the AOI timeline.

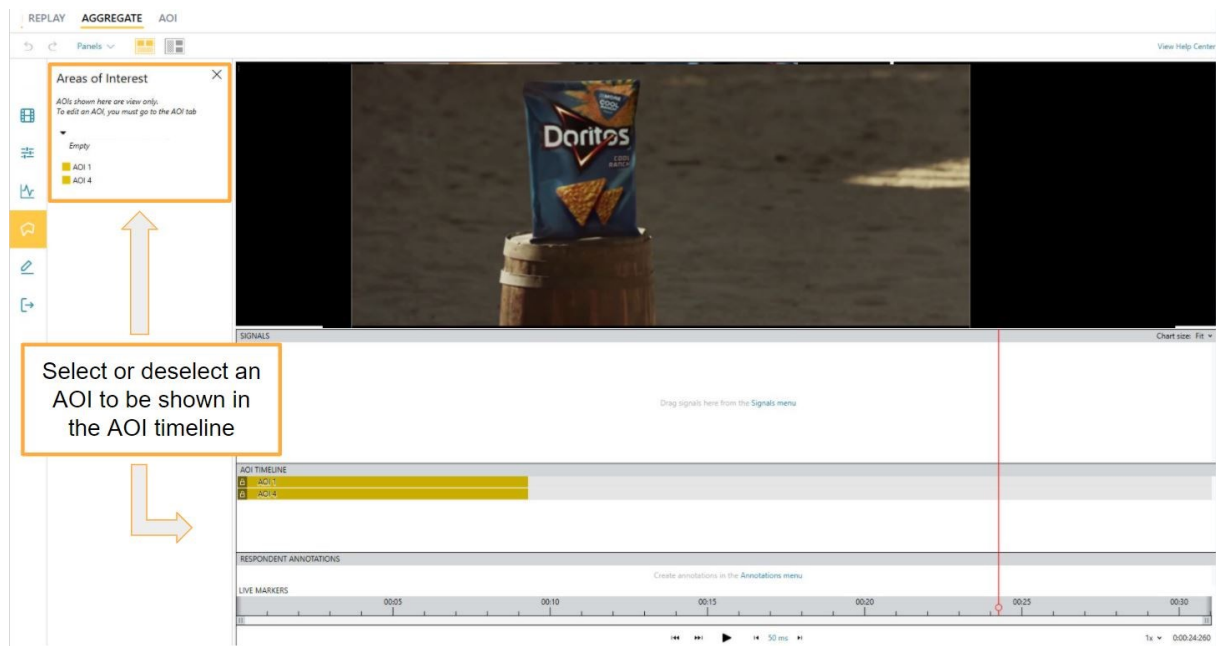

Stimulus annotations

'Stimulus annotations' allow you to [annotate](#) your data from within the Aggregate replay editor. This means that the annotations immediately apply to all respondents in the analysis. However, note that it typically only makes sense to use Stimulus annotation if the annotated events have the exact same time in all respondents.

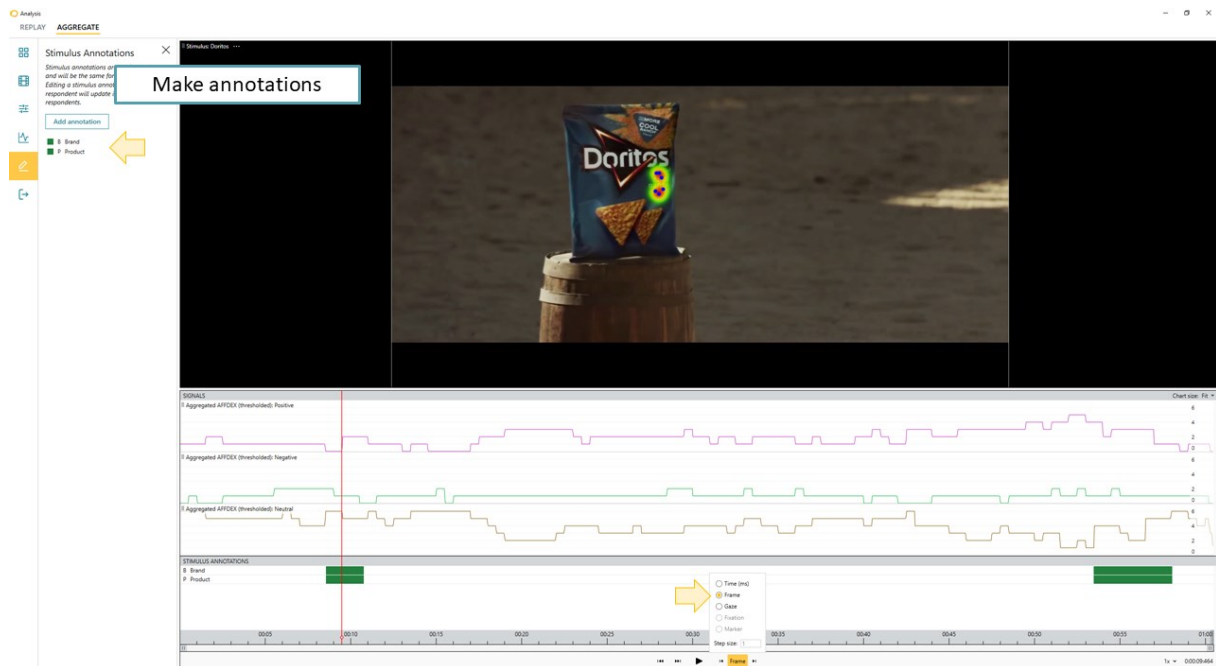

How to add annotations


Add annotations:

1. Navigate to the Stimulus annotations menu
2. Define an annotation by clicking 'Add annotations', and start annotating

Unlike in the Replay editor, in the Aggregate replay editor it is not possible to visualize markers.

Export

The 'Export' menu is finally also very similar to the Export menu in the Replay editor, however, with the additional possibility to export Heatmaps.

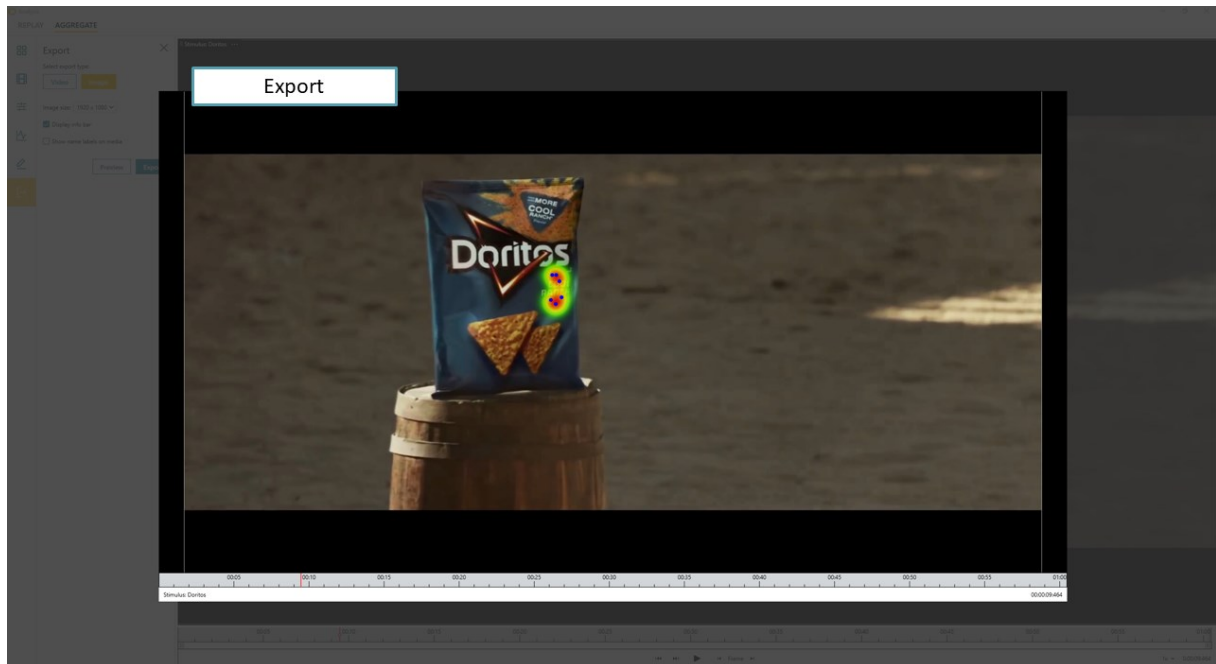

Preview of a heatmap export

Export an image (.png) or video (.wmv, .mp4):

1. Navigate to the Export menu
2. Select Video, Image or Heatmap
3. Click Export

Good to know (features)

A few more features are good to be aware of when operating the Aggregate replay editor. However, these are essentially the same features as in the Replay editor.

See also: [Replay\(Help center article\)](#)

For a comprehensive overview on (1) how to play/pause the data at regular speed, or in fast/slow motion, (2) how to navigate through your data with precision (3) how the zoom, and (4) how to mute the audio, we refer once again to the Help Center article about the Replay editor (see also: [Replay](#)).

