

Table of Contents

iMotions Lab > iMotions Exports

Export: Data Quality Statistics	2
---------------------------------------	---

Export: Data Quality Statistics

This help center article explains how to export data quality statistics.

What is the data quality statistics export about?

The data quality statistics export saves a file (.txt) with basic statistics about the collected data. That includes the *effective* sampling rate of the devices used, and the *percentage* of data obtained; the latter is the inverse of the *percentage* of data missing.


	A	B	C	D	E	F	G	H
1	#Study Name : ...							
2	#Study Date : ...							
3	Respondent_Name	Stim_Name	Shimmer ECG_SamplingRate	Shimmer ECG_Quality	Shimmer GSR_SamplingRate	Shimmer GSR_Quality	Eye tracking_SamplingRate	Eye tracking_Quality
4	R1	Welcome	485	-1	122	93	61	100
5	R1	Instructions	508	-1	127	99	60	100
6	R1	Slide1	499	-1	125	98	60	100
7	R1	Dilemma1	508	-1	127	98	60	90
8	R1	Slide2	505	-1	126	98	61	100
9	R1	Dilemma2	494	-1	124	96	60	100
10	R1	Slide3	504	-1	126	98	60	100
11	R1	Dilemma3	504	-1	127	95	61	100
12	R1	Thank you!	498	-1	125	93	61	94
13	R2	Welcome	487	-1	122	92	60	100
14	R2	Instructions	500	-1	126	95	61	100
15	R2	Slide1	482	-1	121	98	61	100
16	R2	Dilemma1	492	-1	123	96	61	92
17	R2	Slide2	501	-1	126	98	61	100
18	R2	Dilemma2	504	-1	126	96	60	100
19	R2	Slide3	498	-1	125	98	61	100
20	R2	Dilemma3	505	-1	127	95	61	93
21	R2	Thank you!	491	-1	123	93	60	98
22								

These data quality statistics are also visualized within the iMotions software itself.



How to export it?

The data quality statistics export is available *immediately after a recording has ended*, and is exported from the page 'Recorded data'.

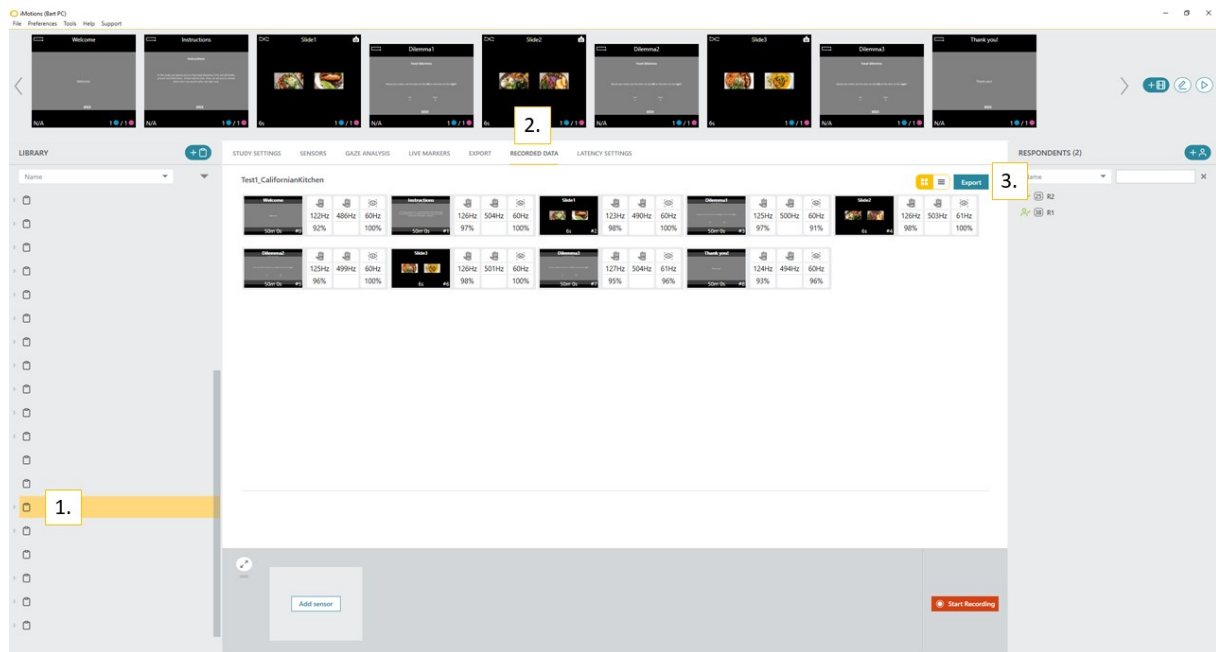

In order to do that (yellow path)

1. Click the name of the study in the Library
2. Navigate to the tab 'Recorded data'
3. Click 'Export', and save the .txt file in the desired folder.

If you now explore the folder, you should find a .txt file with the ending '_DataQuality.txt'.

