

Table of Contents

iMotions Lab > Study Setup > Study Settings and Options

One-Screen Calibration: Configurations and Settings	2
---	---

One-Screen Calibration: Configurations and Settings

This provides a single screen to review sensor connection and participant position.

In iMotions 5.5 and later, the One-Screen Calibration (Compact Calibration) provides the status of all relevant sensors on one single screen, while blending out user interface elements that are irrelevant for determining if all sensors are connected and respondents are positioned correctly.

In some scenarios this functionality allows respondents to “self-calibrate” with limited supervision from an operator or terminate a study if calibration / sensor criteria are not met.

Please note that no live view is available during one-screen calibration in order to minimize any feedback respondents are exposed to. Further, troubleshooting most likely requires the operator to terminate the recording session and cancel the calibration screen. In this case, operators still can select standard calibration to identify and solve any issues.

1. Prerequisites

One-screen calibration has been optimized for single screen setups. It shouldn't be used in dual screen scenarios.

2. Select Calibration Mode for Studies

You can select the calibration mode at the study level. In this case, all respondents (unless configured differently) will undergo the calibration mode set on the study level:

“Regular” calibration mode (the default mode, also set for studies collected in iMotions 5.4 and earlier) triggers the same behavior as previous versions of iMotions, while “Compact” calibration results in the new “One Screen Calibration” mode.

This also affects respondents which are imported using a test plan (since this parameter is not available as a test plan respondent property).

LIBRARY

STUDY SETTINGS SENSORS EXPOSURE STATISTICS RESPONDENT STATISTICS

Date

Study 14-12-15 11h23m

Study 14-12-15 11h23m-1

Study 11-12-15 15h16m

New study

[GSR base] ET + GSR + ECG + FACET

Trawell Advertisements

Name Study 11-12-15 15h16m

Resolution 1920x1080 Custom Resolution

Stimuli Order

☐ Randomize stimuli except for the locked stimuli

☐ Order stimuli as displayed in the stimuli navigator

☒ Load stimuli order from a test plan

<Browse for Test plan> Browse...

Description

Calibration Mode ☐ Regular ☒ Compact

3. Select Calibration Mode for Respondents

For each newly created respondent users can select either “Regular” or “Compact” calibration mode. The “Compact” calibration mode will trigger the “One Screen Calibration” mode.

RESPONDENT DETAILS RECORDED DATA

RESPONDENTS (1)

Name

John

Gender ☒ Male ☐ Female Age 0

Group Default

Comments

Calibration Mode ☐ Regular ☒ Compact

Update Undo

@ John

4. One-Screen Calibration Window


1. Select the "Compact" calibration mode at the study or respondent level.
2. Click the "Record" button to start the study.



3. If eye tracking has been enabled, the respondent will undergo the eye tracker calibration.
4. Upon successful calibration, the One-Screen Calibration window is displayed, visualizing the quality of all enabled sensors (e.g., FACET, Shimmer GSR, eye tracking). The color bar below each sensor indicates signal quality ("green" = optimal; "red" = non-optimal).

Mike

Edit Details

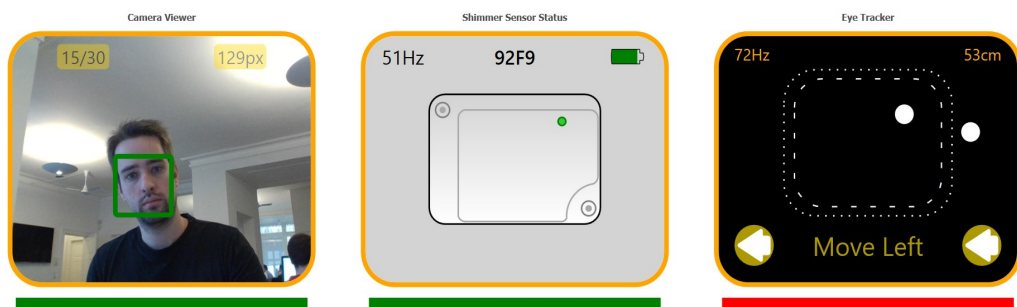

5.

Back

Override

✓
6. In the top left corner the name of the respondent is visible. You can also change the details by pressing "Edit Details"

Mike

Hide Details

Name: Mike


Age: 28

Gender: ☒ Male ☐ Female

Group: Default

Save

Cancel

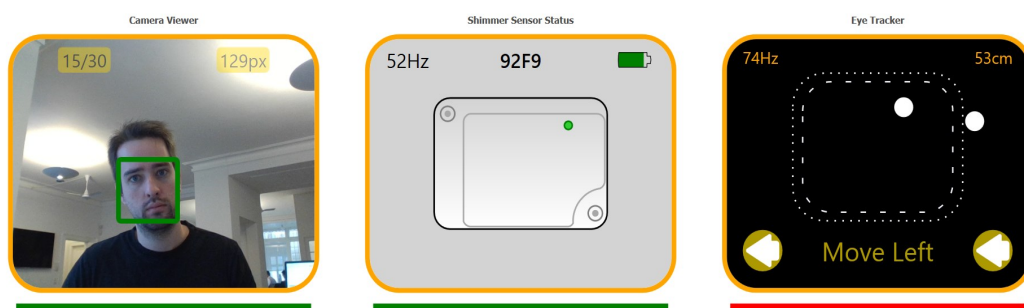


7.

Back

Override

Continue
8. In "Edit Details" you will see "Name" "Age" "Gender" "Group" which can be

changed directly here in "One screen calibration" by pressing "Save". Changes made will be saved to the study in iMotions UI. To undo any changes to the respondent press "Cancel".

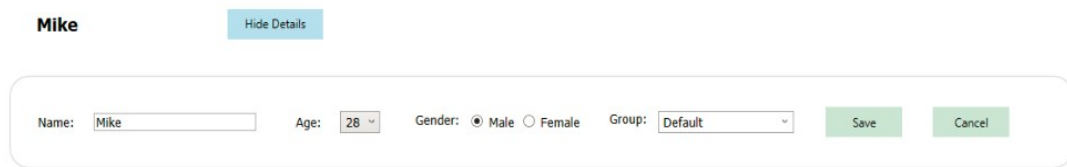

The screenshot shows a calibration form for a respondent named "Mike". At the top left, the name "Mike" is displayed next to a "Hide Details" button. Below this, a form contains fields for "Name" (filled with "Mike"), "Age" (a dropdown menu showing "28"), "Gender" (radio buttons for "Male" and "Female", with "Male" selected), and "Group" (a dropdown menu showing "Default"). To the right of these fields are two green buttons: "Save" and "Cancel".

9. If all sensors are "green", you can click "Continue" to proceed with the data collection. If one or more sensors are "red", you can click "Back" to repeat any sensor setup and calibration or "Override" to proceed (not recommended).



Auto advance

In Global Settings you can set up the One Screen Calibration to autoadvance to the study if eye tracking passes a quality threshold.

Global Settings

Eye tracking Calibration Slide Show Stimuli Results Video Sensors API

☐ Eye Tracking Server

Model: Tobii 1750 ☒ Forward samples

☒ Remove low quality data

Test Connection

Eye Tracker Calibration

Number of Calibration Points: 9

Points Order and Location: Use Default Customize

Calibration point Size: Medium

Calibration point color: [Color Picker]

Calibration background color: [Color Picker]

Calibration point speed: Medium

Calibration point Visualization: ☒ Default ☐ Custom

Image/Video: [Image/Video Picker] Add Clear

☒ Calibration points transition

☒ Continue on poor calibration

☐ Use standard calibration for scene recording


☒ Start slide show after calibration with minimum quality: Excellent

OK Cancel

Override


If the respondent does not pass the threshold set in Global Settings for Calibration the operator can chose to override and continue to the study. This will typically happen in case the calibration has failed a number of times but the respondent's data is still deemed useful regardless of the eye tracking quality (e.g. if data from other sensors can be used). Here is how the screen looks like when respondent did not pass the calibration:

The Gaze Calibration Quality Is Below The Selected Threshold

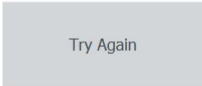

Try Again

The operator can use the key combination for Next-Slide. If you do not know the key combination for Next-Slide you can find it under Global Settings->Slide Show. The default value is SHIFT+SPACE.

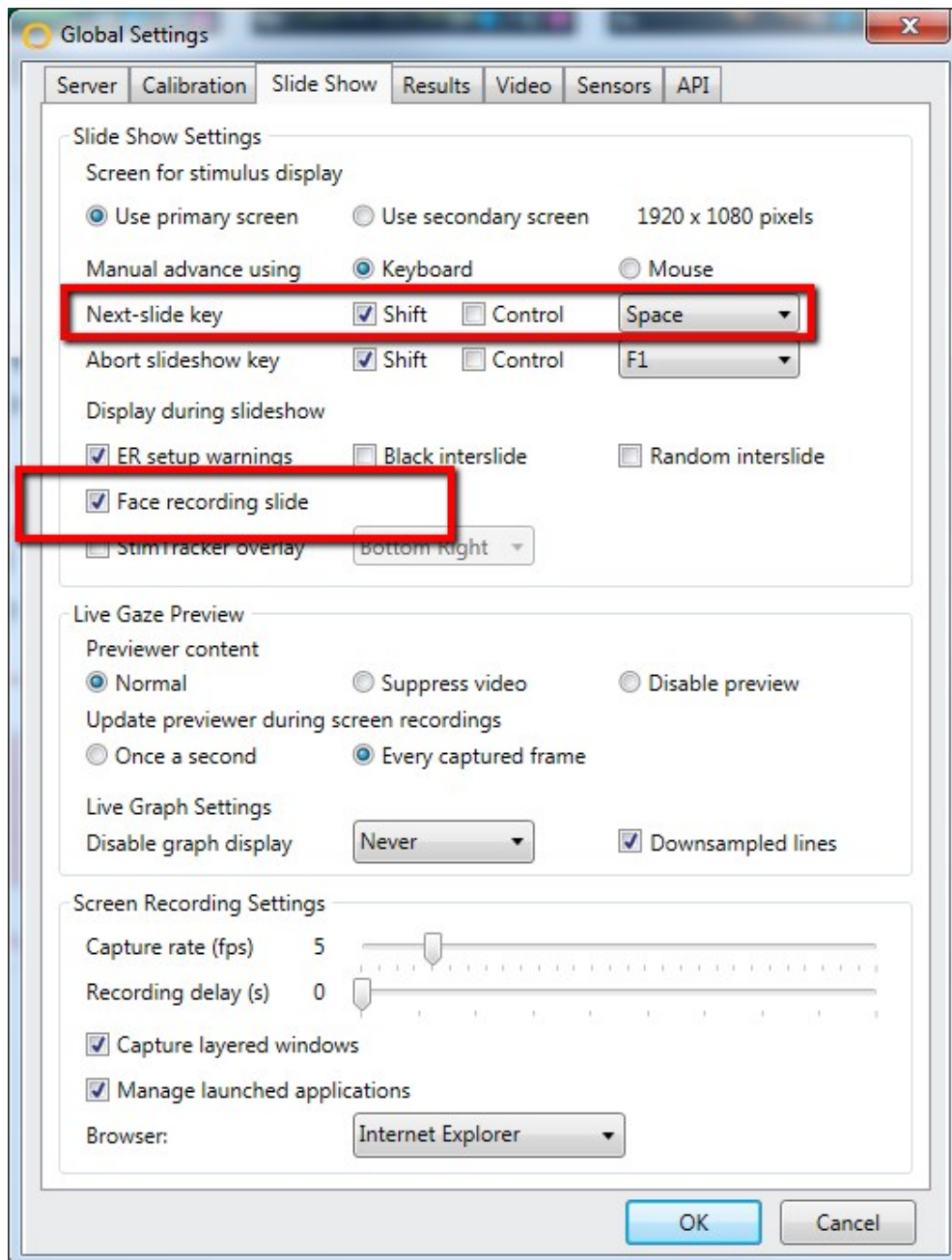

The face recording slide or placeholder slide will be displayed in single screen mode if so requested in global settings, this would help you if you are running a single screen setup and wish to hide the operator screen.

