

Table of Contents

iMotions Lab > Analysis > Signal Processing (R Notebooks) > Facial Expression Analysis

Re-Encoding of Incompatible Video Stimuli	2
---	---

Re-Encoding of Incompatible Video Stimuli


This article is valid for iMotions version 7.0 and onwards.

Overview:

This article describes how you can easily re-encode video stimuli that are not working properly in iMotions, e.g. lagging, flickering, or with poor audio. Note, this feature will not convert video files to different formats, but will re-encode your video stimuli to iMotions' default and replace it with the original stimulus.

How to identify erroneous video stimuli:


Ideally, erroneous video stimuli should be identified before the actual data collection step. The easiest way to make sure that your video stimuli are running properly in iMotions, is to use the preview function. After adding stimuli in iMotions, your study can be previewed by clicking the preview button as illustrated below.

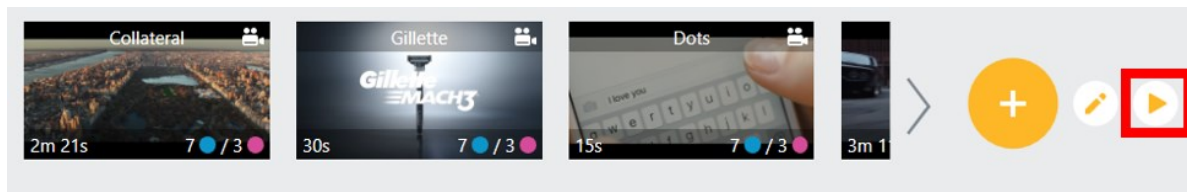

Note, re-encoding of your video stimuli can also be done after data has been collected in iMotions.

How to re-encode video stimuli:


To re-encode your video stimuli, right-click on the study containing the erroneous stimulus and press 'Convert videos for compatibility' in the occurring window as illustrated below. This may take a moment and will only work for video stimuli.

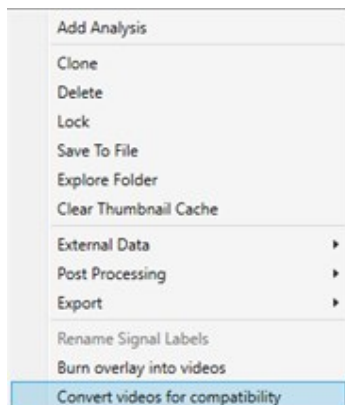

Advanced feature use for face recordings:

In case you recorded videos outside of iMotions and experience erroneous videos when you import them for facial coding with the [AFFDEX post-processing](#) feature, the following can be done:

1. Create a study in iMotions and [add the recordings as video stimuli](#).
2. Re-encode the recordings with the 'Convert videos for compatibility' functions as described above.
3. Locate the re-encoded videos from the study folder by right-clicking on the study and selecting 'Explore Folder' as illustrated below. The converted videos will be located in the 'Stimuli' folder.

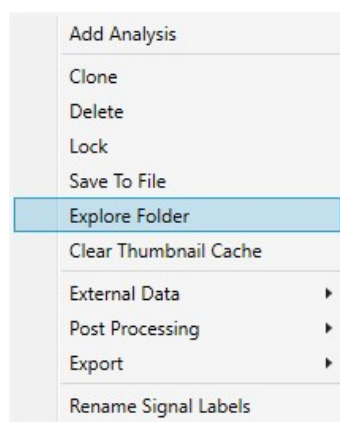

4. Set up your study for [AFFDEX post-processing](#) using the re-encoded videos.

