

# Table of Contents

iMotions Lab > Analysis

|                                   |   |
|-----------------------------------|---|
| Replay and Analysis Windows ..... | 2 |
|-----------------------------------|---|

# Replay and Analysis Windows

## How to Review and Analyze Data in iMotions Replay Windows


iMotions Replay Windows allow you to review and analyze your recordings. Each window displays the timeline of your recording and a playback navigator, along with displayed panels and menus. Each Replay window contains different menus, panels, and functions (described here).

- The Respondent Replay is available immediately after the recording is completed. This displays sensor data, stimuli, and camera recordings.
- Additional Analysis Replay windows are available, once an study analysis is created.
  - This Replay window displays the same content as a Respondent Replay but also includes analysis results.
  - The Aggregate window displays group results, if they apply to the stimulus (e.g. video/image media files, Scenes).
  - The AOI window is used to create AOIs and calculate AOI Metrics (only available if your study includes eye tracking data).


## Respondent Replay Window

### Opening Respondent Replays

1. Select Study that contains Respondent Recording(s)
  2. Select Respondent
  3. Click Play button for any stimulus to open the Respondent Replay
  4. Use the Menu Icons and Panels options to select content to display (empty panels also provide buttons that will open the corresponding menu)
-

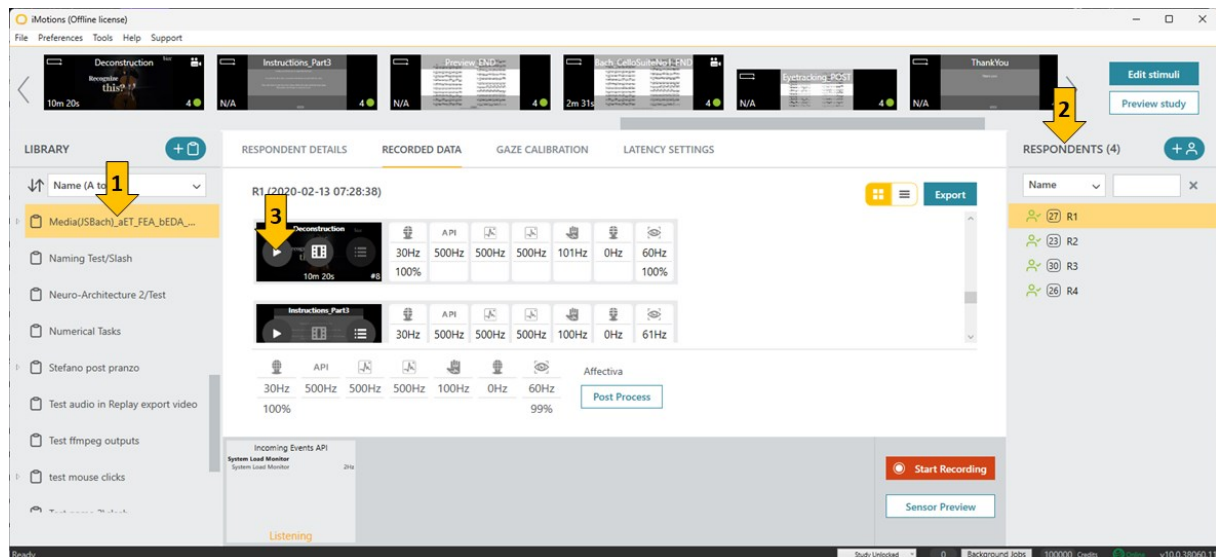

Opening the Respondent Replay Window

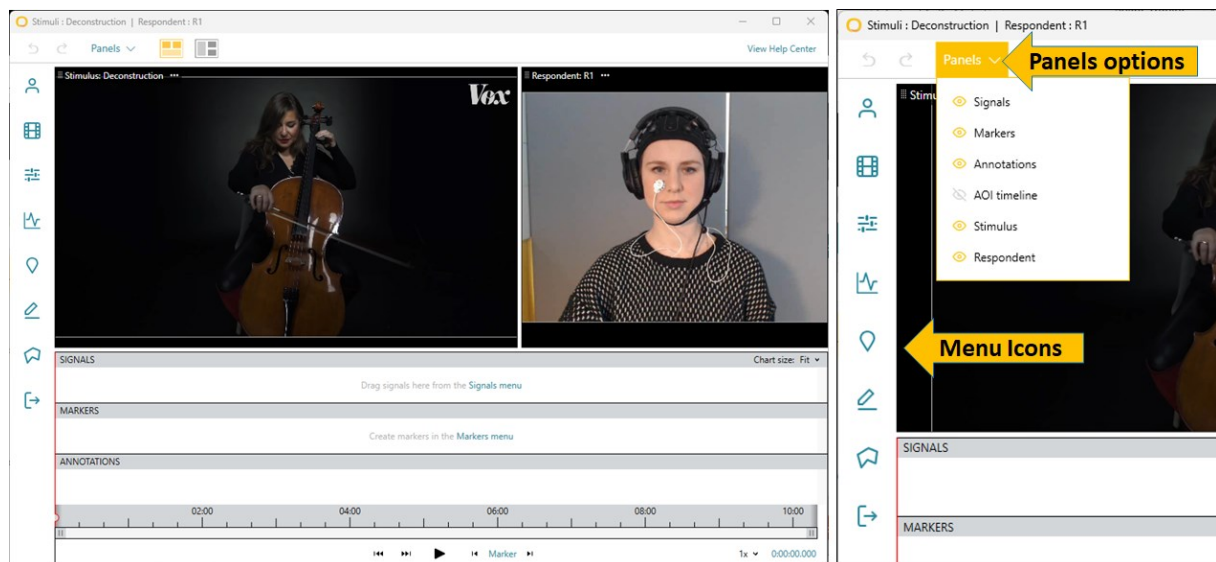

Respondent Replay (left) with Menu Icons and Panel options (right)

## Respondent Replay Menus

1. Respondents - List of recorded respondents
2. Recordings - List of stimuli associated with recordings
3. Settings - Display settings for applicable overlays, allowing visualization of Gaze (from eye tracking data), Mouse (showing cursor location), AOIs, or Annotation comments.
4. Signals - List of sensors that expand to show available data streams. You can select which data will appear by clicking on the data
5. Markers - List of Event Markers
6. Annotations - List of Annotations
7. Areas of Interest - List of Areas of Interest
8. Export - Options for exporting images or videos from your replay




Menu Icons

## Creating an Analysis

Additional Analysis Replay windows are available, once an study analysis is created. General instructions are provided here. Please refer the section of our Help Center about using specific [R Notebooks](#) or [Areas of Interest \(AOIs\)](#) in your analysis.

1. Click the Create Analysis button, which appears when you hover over a study in your iMotions Library.

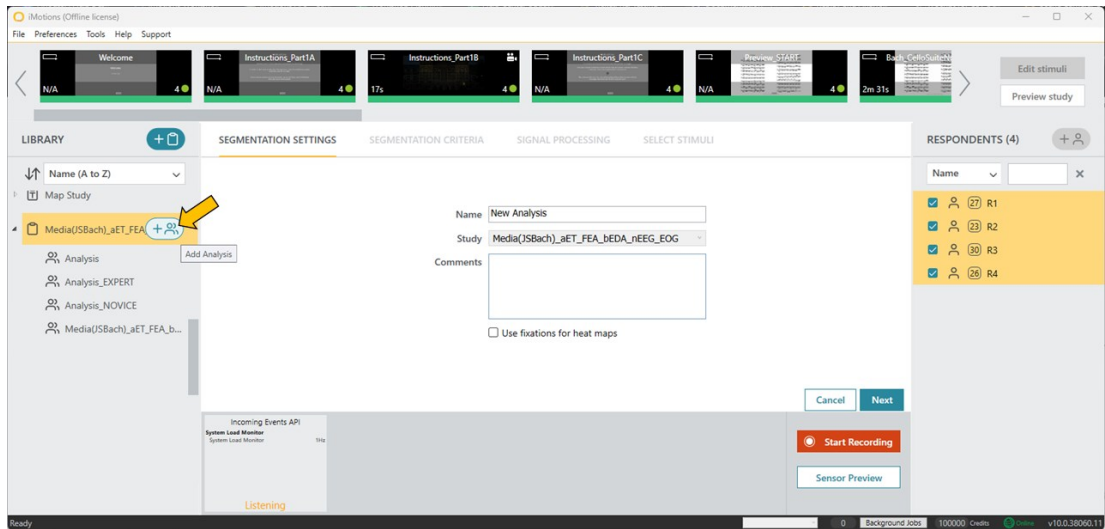

Create Analysis button appears when hovering over a study name


## 2. Enter Analysis settings

- a. Segmentation Settings - Options for Name, Comments, and 'Use fixations for heat maps' setting.
- b. Segmentation Criteria - Options for selecting Respondents based on segmentation criteria for Gender, Age, Flow, Group, and Advanced Segmentation options (e.g. data quality metrics).
- c. **\*\* Signal Processing \*\*** - Workflow options for your analysis. You can create an analysis without any 'Active flows'. Note: [AOI analysis](#) done, separately.
  - i. Devices - List of sensor devices included in your study
  - ii. Algorithm - List of R Notebook workflows for the selected device
  - iii. Flow name - Option to name the signal processing flow (must be unique from other flow names in this analysis)
- d. Select Stimuli - List of stimuli to select and include in the analysis (all stimuli, by default).

Signal Processing settings define your analysis workflow via R Notebook selection and settings

## Analysis Replay Windows

Once an analysis is created/selected, you can click/open the Replay, Aggregate, or AOIs windows. The Gaze Map window allows you to replay and analyze recordings with eye tracking data, which is covered in a designated Help Center article (due to difference in the interface and workflow).

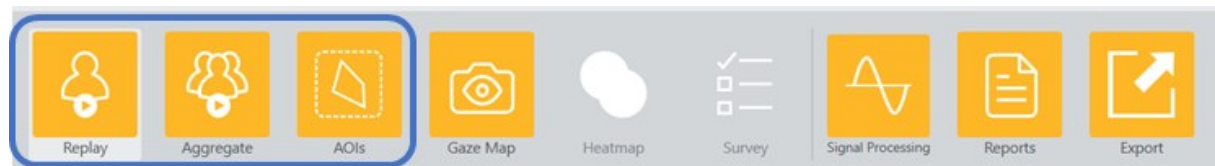


Analysis Replay Options for (analyzed respondent) Replay, Aggregate, and AOIs windows

## Replay

This Replay window (within an analysis) displays the same content as a Respondent Replay but also includes additional information and options in Signals. Results are generated by any Signal Processing that has been performed in the Analysis.

This Signals Menu includes raw sensor data, analyzed data, and summary metrics that are available.


Note: A 'Respondent Replay' opens from the Study level (i.e. the Study is selected vs. an Analysis within the study) and will not contain analyzed data or summary metrics.

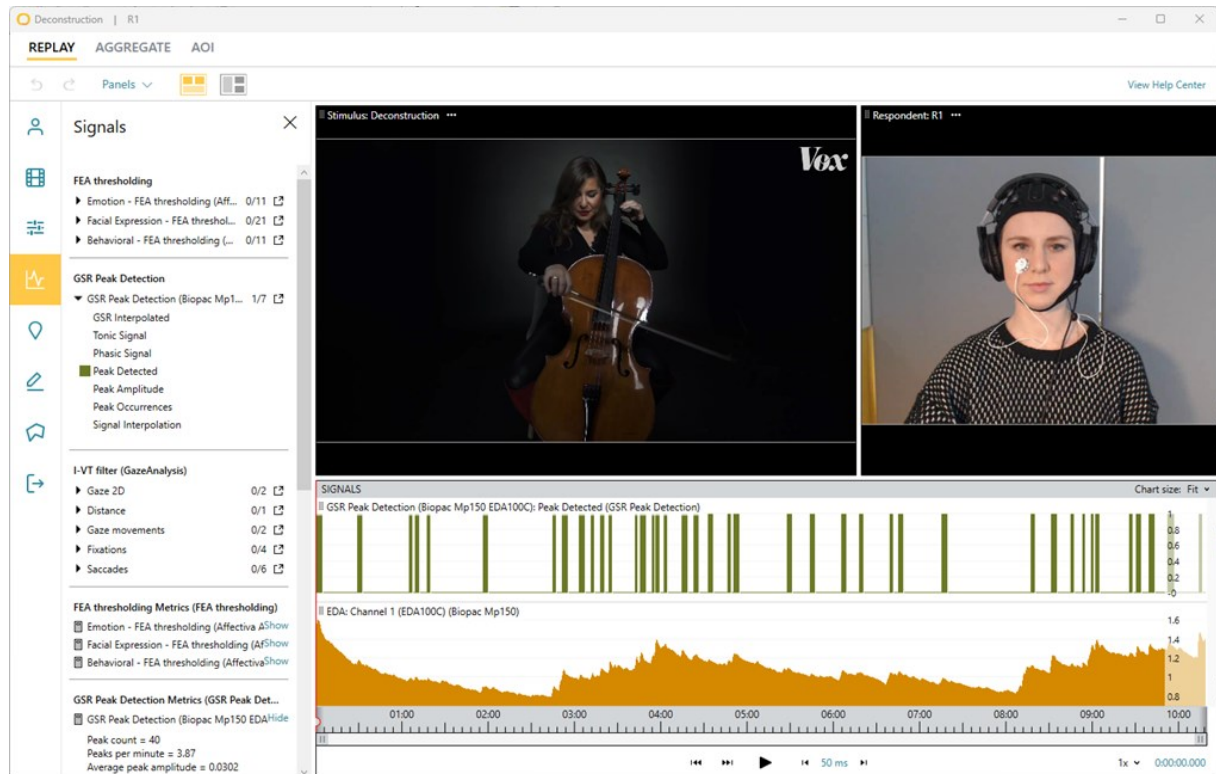


Analysis Replay Signals shows available signals (including analyzed data) and display settings for the Signals Panel

## Aggregate

- The Aggregate Replay displays group results, if they apply to the stimulus (e.g. video/image media files, Scenes). If Aggregate results do not apply, the Aggregate button will be greyed out and unavailable.
- Aggregate Signals can be selected in the Signals Menu and displayed in the Signals Panel.
- Aggregate eye tracking data can be displayed over the stimulus video (in the Settings Menu). These settings control the appearance of gaze data and whether AOIs are displayed. Some of these options are also available by right-clicking in the Stimulus Panel.

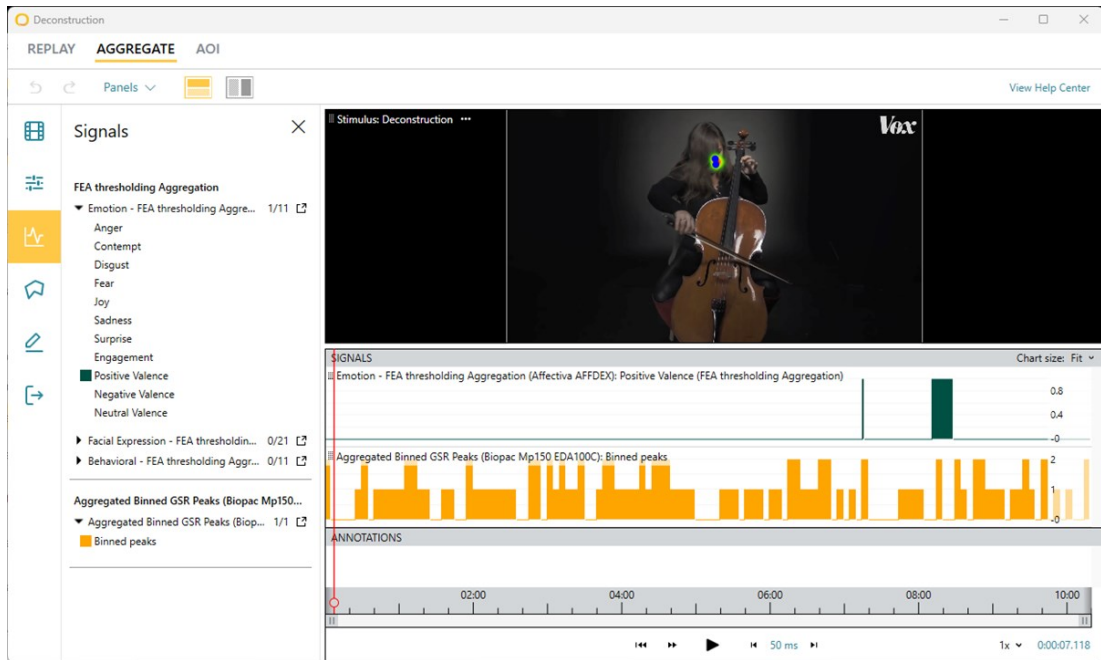

Aggregate Signals Menu shows available signals and display settings for the Signals Panel

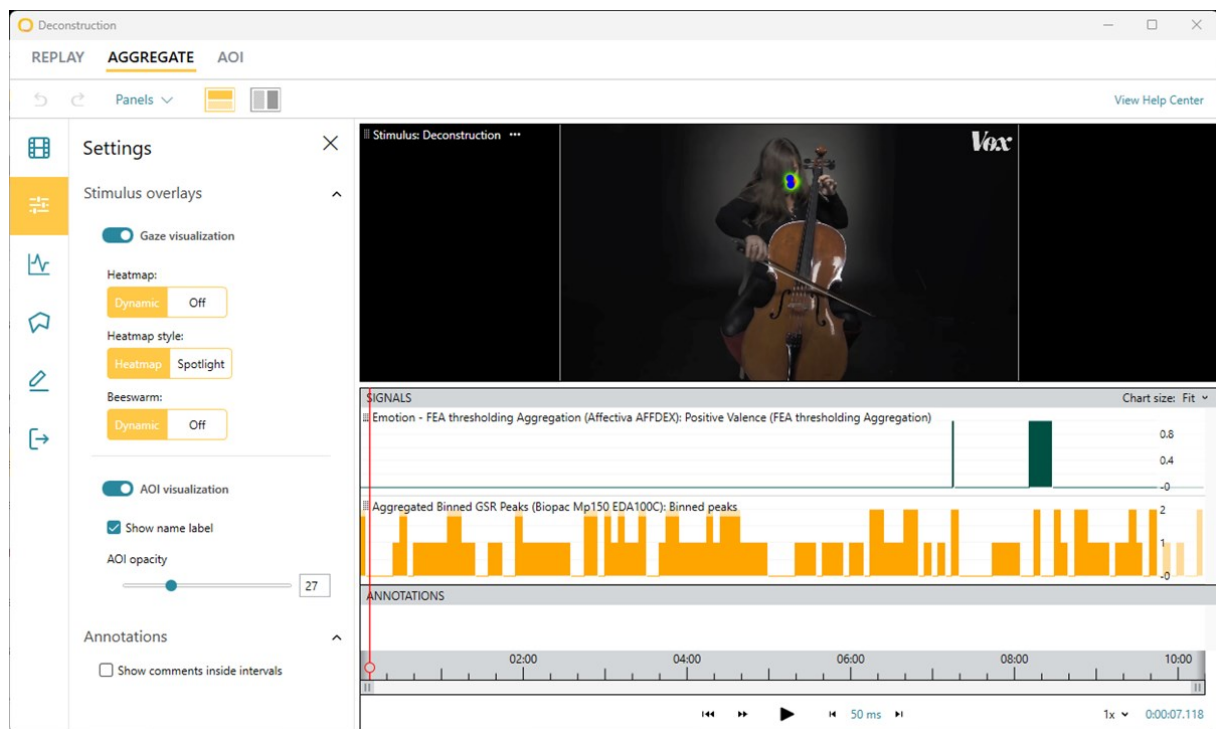


Aggregate Settings Menu showing gaze visualization options


## AOIs

- The AOIs Replay is used to create AOIs and calculate AOI Metrics (only available if your study includes eye tracking data). Detailed information is available in our Help Center article about [Areas of Interest \(AOIs\)](#).
- The 'Areas of Interest' menu is used to create AOIs and edit properties.
- The 'AOI Metrics' menu is used to select the metrics that will appear in the 'AOI Metrics' Panel. Clicking on any AOI metric will toggle whether it is

displayed. Information about each AOI metric in the AOI Metrics Menu, by hovering over each metric then the "?" icon.

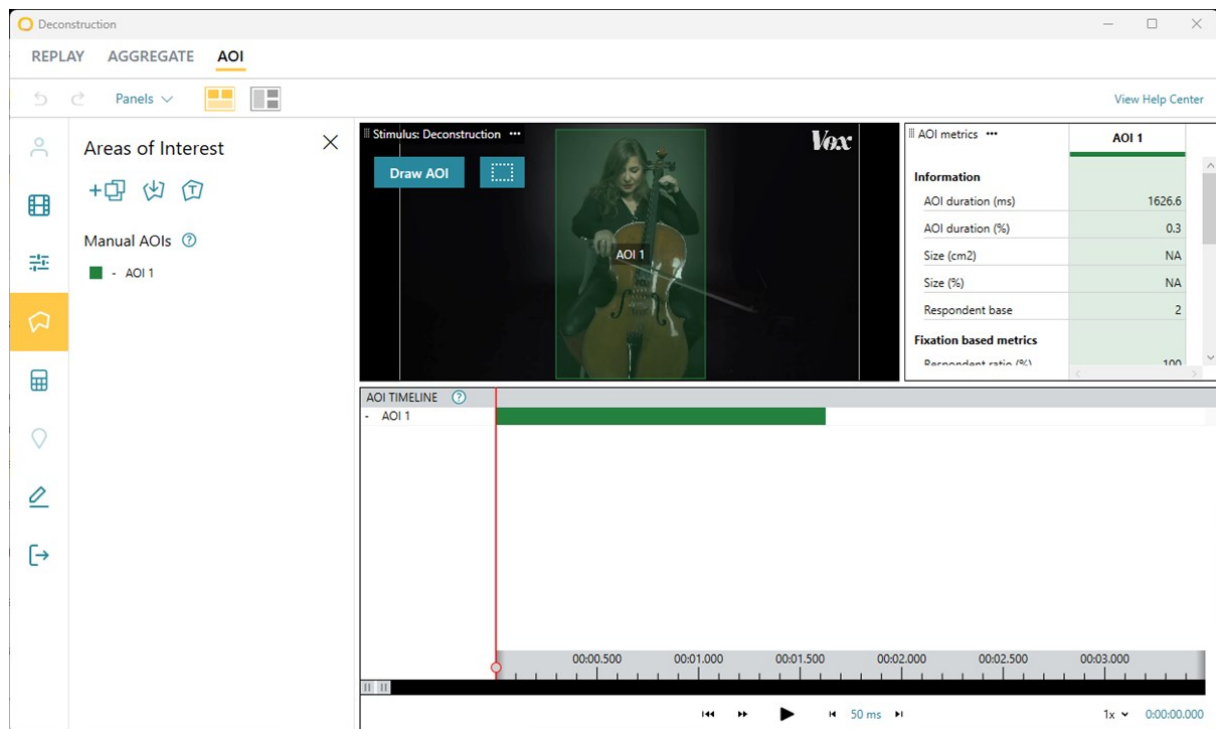


AOIs are created in the AOI replay window using the 'Draw AOI' button or 'Areas of Interest' menu

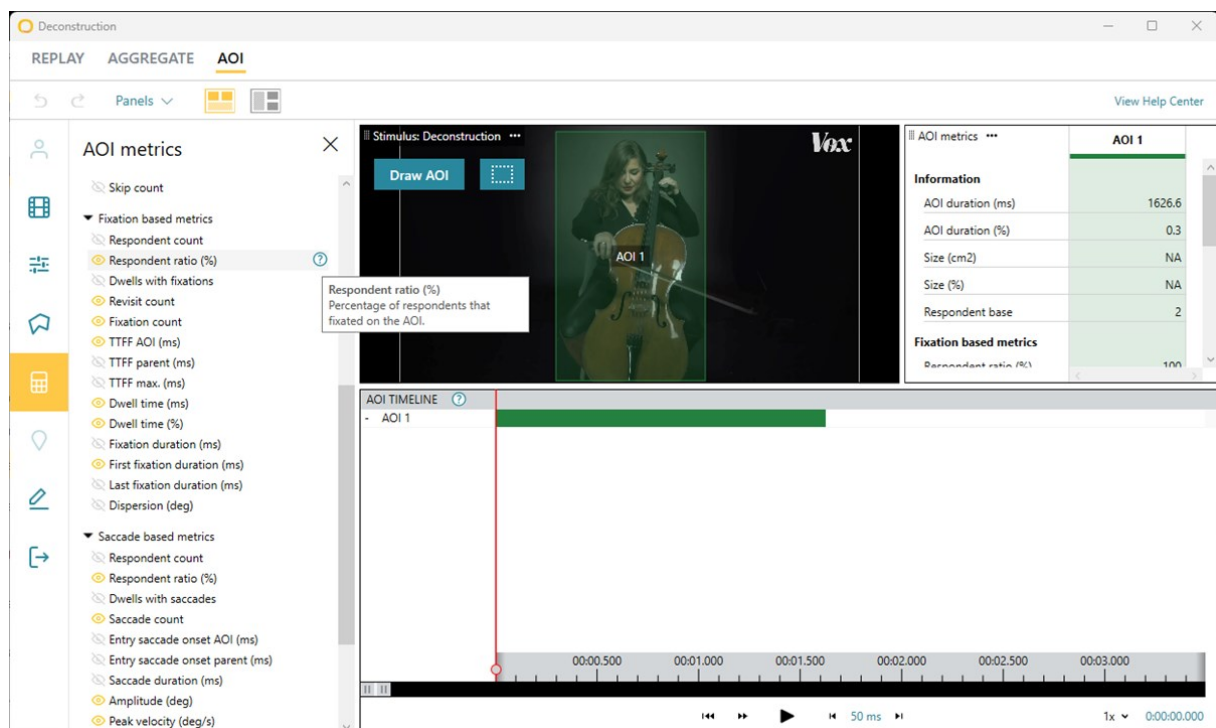

The 'AOI metrics' menu is used for initial calculation of AOI metrics and selecting metrics shown in the 'AOI metrics' panel


# Navigating Replays (Playback functions)

The Recording Timeline and Playback Navigator is available at the bottom of all replay windows. You can use the mouse or keyboard to begin playing the recording or select time points for review and analysis.

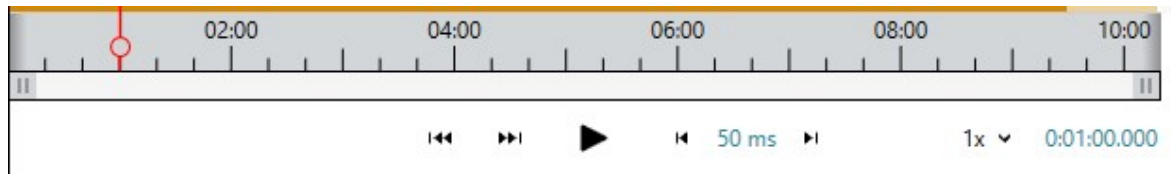

Recording Timeline and Playback Navigator

- The navigator displays (from left to right) the Start/End of Recording buttons, Play button, 'Step' buttons, Playback speed options, Time point.
- The red cursor/line indicates the time point currently shown from the recording. This line extends to any time series panels that are displayed (e.g. Markers, Annotations, and AOI Timeline)
- You can use your mouse to click on any time point or drag this cursor.
- You can also use the keyboard to Play/Pause (with the spacebar) and step through the replay (with the arrow keys).
- Step settings are available between the Forward and Backward Step buttons. This includes can be set to a specific number of milliseconds or frames or to event markers.
- Playback speed can also be adjusted (bottom right). 0.25x, 0.5x, 1x, and 2x speeds are available.
- You can find the current time point (bottom right corner). Clicking here allows you to enter a specific time point and provide an option to display time-of-day (e.g. 12:00 PM) instead of the time in the recording.

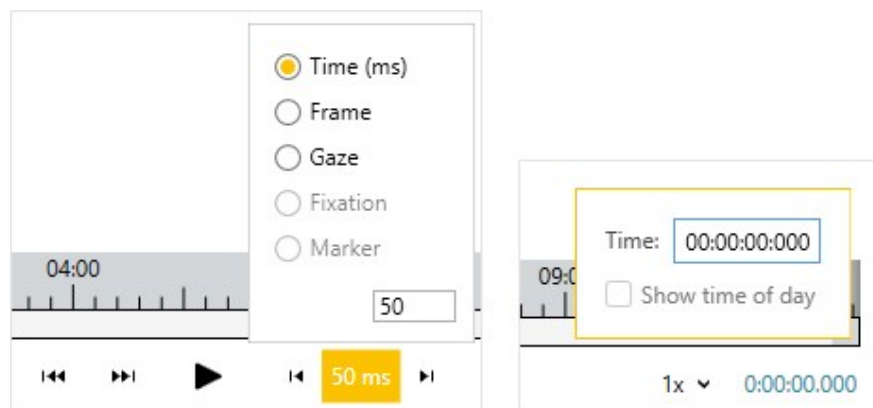

Playback 'Step' options and Time input/display option

