

Table of Contents

iMotions Lab > Sensor Setup > Screen-based Eye Trackers (SBET)

Setting up Participant-guided Eye Tracker Calibration	2
---	---

Setting up Participant-guided Eye Tracker Calibration


This article walks through the setup and use of participant-guided eye tracker calibration.

Participant-guided eye tracker calibration has calibration steps presented on the participant screen (in iMotions 10). At the start of data collection, the eye tracking preview widget will be shown on the participant's screen. Calibration begins after clicking the "Continue" button with the mouse or pressing the Next-slide keyboard buttons (defined in the iMotions settings > Slide Show tab).


1. To enable participant-guided eye tracker calibration, go to Study Settings and select "Respondent guided" under "Calibration mode". This option is only available, if eye tracking is selected in Sensors settings.

The screenshot displays the 'Study settings' window in iMotions. The 'STUDY SETTINGS' tab is active. The 'Study name' field contains 'Respondent guided Calibration'. The 'Description' field is empty. The 'Study type' is 'Screen-based', with a sub-description: 'Stimuli shown on a screen, (i.e. images, videos, websites, surveys, etc)'. The 'Screen' section shows 'Resolution: 2560 x 1080', 'Width: 67.3 cm', 'Height: 28.4 cm', and 'Distance: Auto'. The 'Calibration mode' section has two options: 'Instructor controlled (recommended)' (unselected) and 'Respondent guided' (selected). At the bottom right, there are 'Undo' and 'Update' buttons.

2. At the start of data collection, after pressing the "Start Recording" button, the sensor preview widget(s) will appear on the participant's screen.

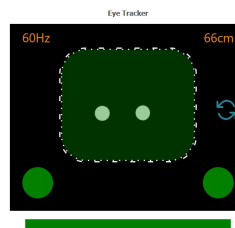


3. When sensors are ready (as indicated by green bars in the widgets), you can click "Continue" to proceed with the data collection. You can also click "Back" to cancel the recording and repeat any sensor setup (or "Override" if a calibration has been done but did not meet the quality threshold defined in Global Settings > Eye tracking)

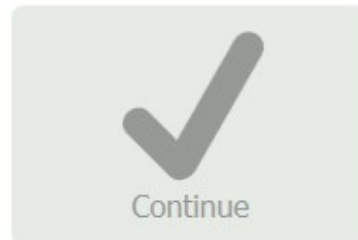

4. Eye tracker calibration will proceed. If calibration quality is not sufficient, the participant is prompted to try again (returning them to the initial screen)

The Gaze Calibration Quality Is Below The Selected Threshold

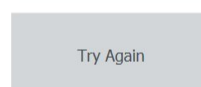

Setting Auto advance threshold

In Global Settings you can set up the One Screen Calibration to auto-advance to the study if eye tracking passes a quality threshold.

The screenshot shows the 'Global Settings' window with the 'Calibration' tab selected. The 'Eye Tracking Server' section is at the top, followed by the 'Eye Tracker Calibration' section. In the 'Eye Tracker Calibration' section, the 'Start slide show after calibration with minimum quality' checkbox is checked and highlighted with a red box. The dropdown menu next to it is set to 'Excellent'. Other settings visible include 'Number of Calibration Points' set to 9, 'Points Order and Location' set to 'Use Default', 'Calibration point Size' set to 'Medium', 'Calibration point color' set to white, 'Calibration background color' set to dark gray, 'Calibration point speed' set to 'Medium', 'Calibration point Visualization' set to 'Default', 'Image/Video' set to an empty dropdown, 'Calibration points transition' checked, 'Continue on poor calibration' checked, and 'Use standard calibration for scene recording' unchecked.

Global Settings

Eye tracking Calibration Slide Show Stimuli Results Video Sensors API

☐ Eye Tracking Server

Model: Tobii 1750 ☒ Forward samples

☒ Remove low quality data

Test Connection

Eye Tracker Calibration

Number of Calibration Points: 9

Points Order and Location: Use Default Customize

Calibration point Size: Medium

Calibration point color: [White]

Calibration background color: [Dark Gray]

Calibration point speed: Medium

Calibration point Visualization: ☒ Default ☐ Custom

Image/Video: [Empty] Add Clear

☒ Calibration points transition

☒ Continue on poor calibration

☐ Use standard calibration for scene recording

☒ Start slide show after calibration with minimum quality: Excellent

OK Cancel

