

Table of Contents

iMotions Lab > Study Setup > Study Stimuli > iMotions Surveys


iMotions Survey Stimuli and Elements	2
--	---

iMotions Survey Stimuli and Elements

This article will walk you through the survey slide designer interface updated to make creating professional looking surveys and questionnaires for your iMotions studies easier.

Note

Please note that survey elements are only displayed properly if Windows display settings are set correctly. In detail, please make sure to set the "Display Scaling" (size of text, apps, and other items) on all screens to 100%. If the display scaling is set to other values, survey elements might not be displayed correctly. Please note that changing the display scaling does not affect how surveys are displayed when running an RDC study

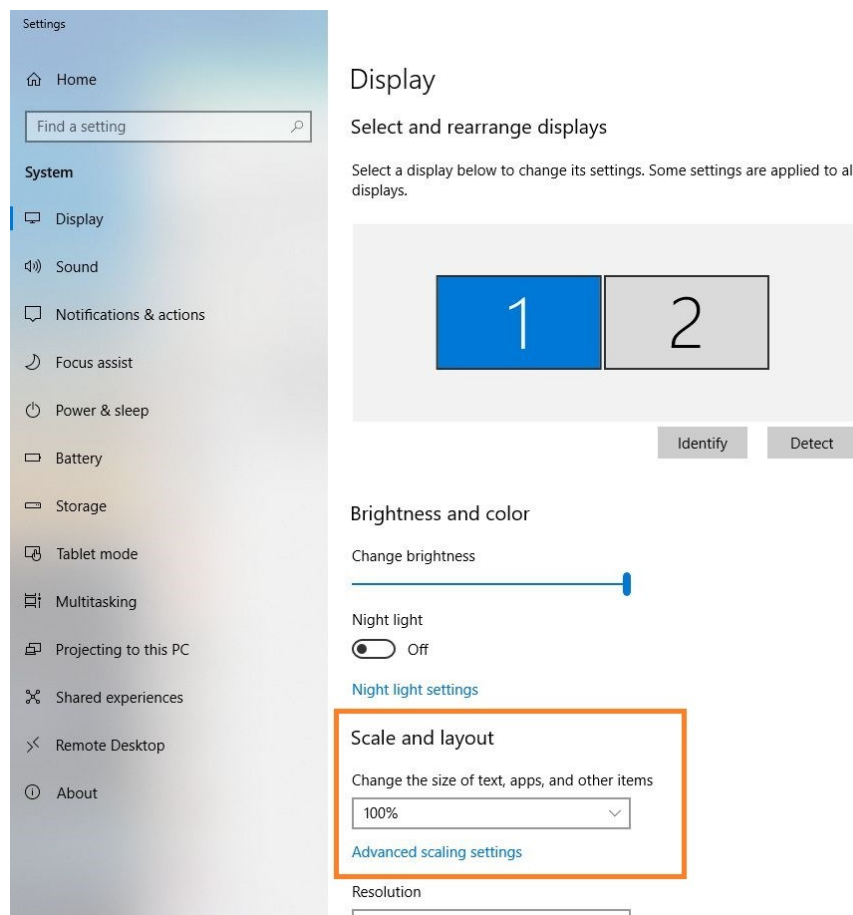

Contents

1. Launching the Survey Slide Designer

2. Designer Interface Overview

1. Work space

2. Elements, Settings, and Templates
3. Properties
4. Preview and Save Bar

3. Element Property Details


1. Rich Text Message
2. Plain Text Message
3. Image Element
4. Scale
5. Single / Multiple Choice
6. Vertical Selection
7. Text Response
8. Key Button Response

4. Survey Settings

5. Templates


1. Creating a Survey Slide

There are two ways to launch the survey slide designer:

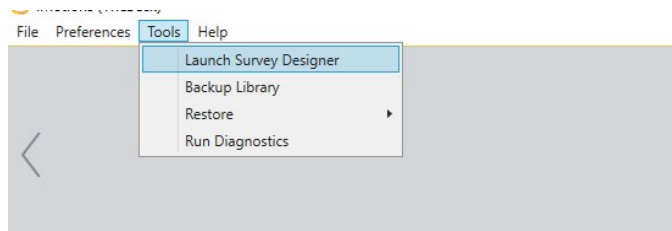


OR

By clicking on Edit Stimuli and adding a Survey

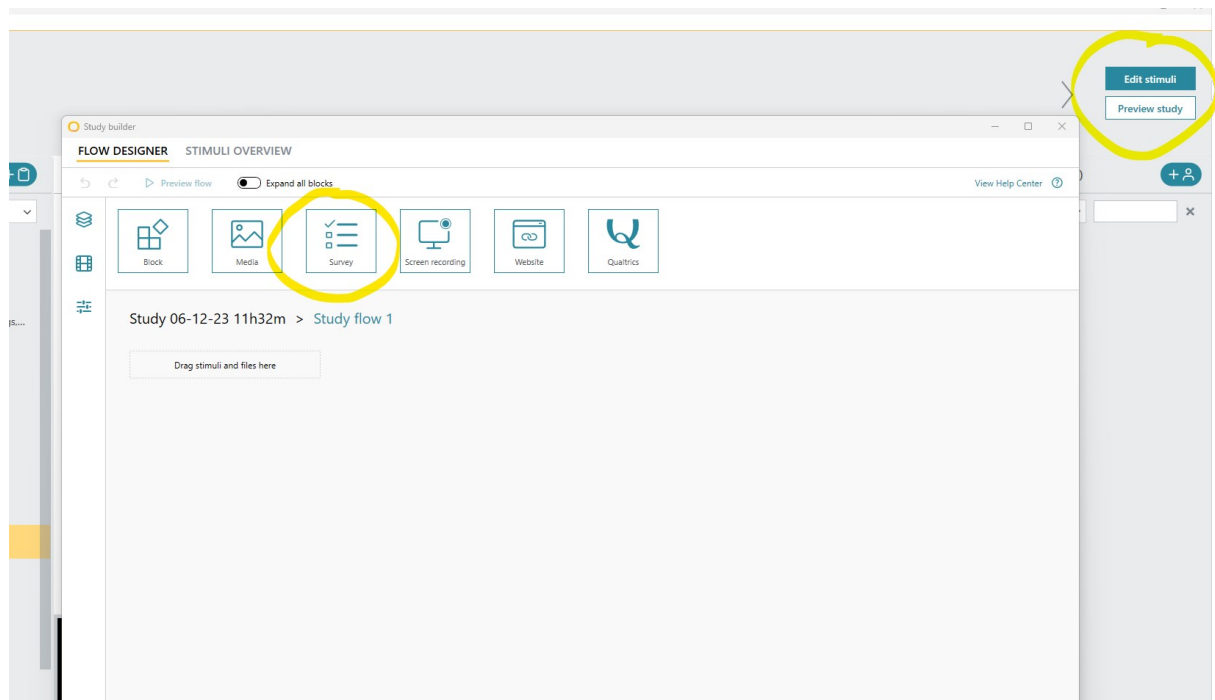


2. Designer Interface Overview

The survey slide designer interface consists of 4 main regions; the work space in the center, the elements, settings, and templates panel on the left, the properties panel on the right, and the preview and save bar at the bottom.

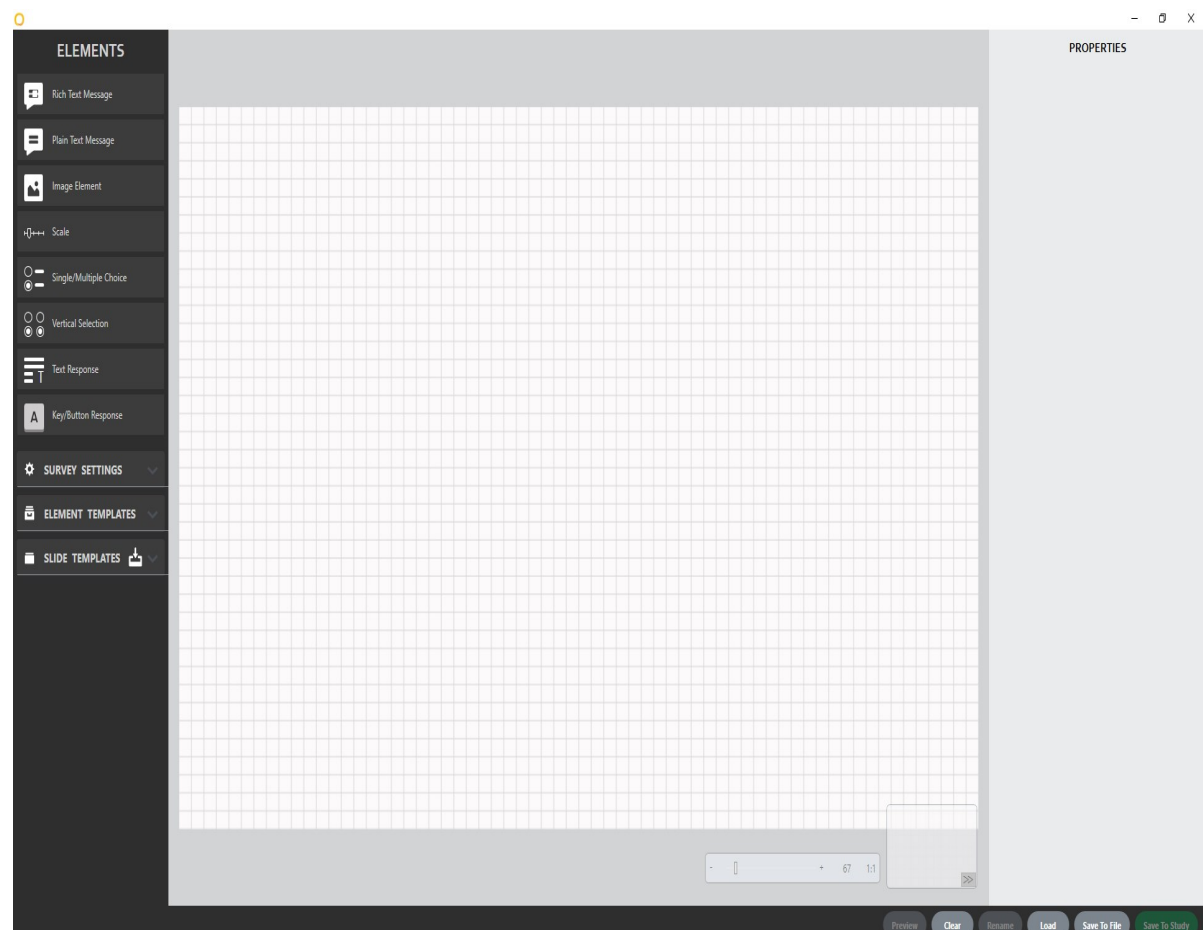


2.1 Work space

The work space is your blank canvas to build your survey on. You can zoom in or out using the tool in the bottom right corner of the work space, or hide the tool as shown.

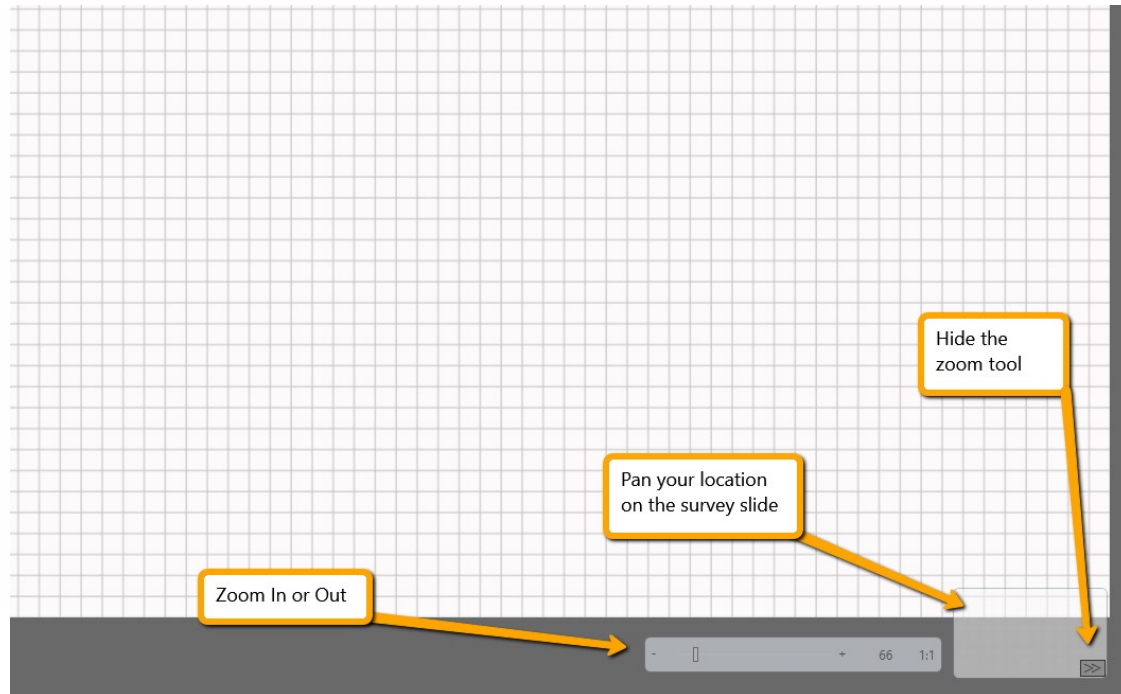

2.2 Elements, Settings, and Templates

The left side bar consists of your survey Elements, settings, element templates, and survey templates. Elements are items you can add to your work space by double clicking on the desired type. Once added, they can be dragged around the work space, customized, deleted, or saved as element templates. Survey settings adjust your work space and overall survey look. Templates make quickly reproducing previously created elements or full survey slides easier.

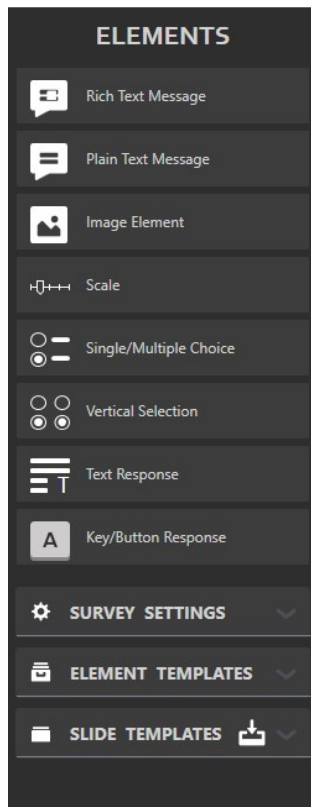


2.3 Properties

Once you have added an element to the survey, selecting your element will bring up properties on the right.

All elements have the following properties:

Center on Slide - Places the element in the middle (horizontally not vertically) of the slide.

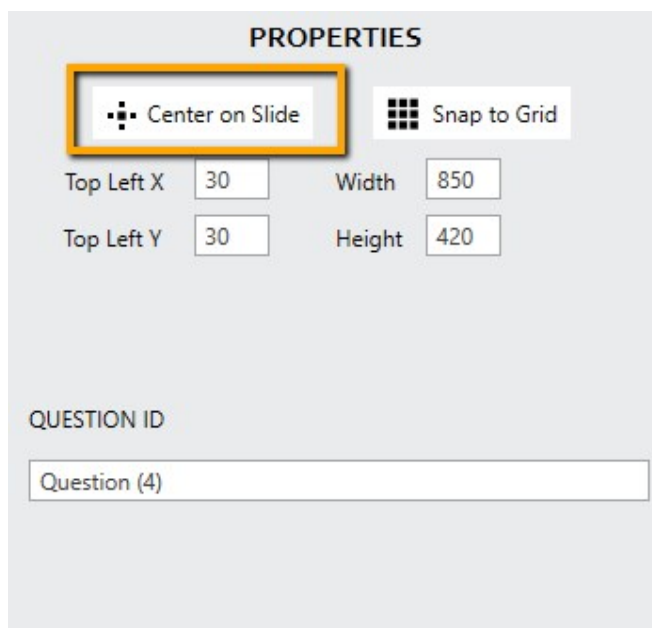

Snap to Grid - aligns the edges of the element with the nearest grid lines in the workspace

PROPERTIES

☒ Center on Slide ☒ Snap to Grid

Top Left X Width

Top Left Y Height

QUESTION ID

Top Left X & Y - set the location of the top left corner of the element

PROPERTIES

☒ Center on Slide ☒ Snap to Grid

QUESTION ID

Width & Height - adjust the size of the element

Question ID- A chosen title for the question. Note that this is what will appear in your survey results in association with the responses so make sure you will recognize your question in the export by this name.


Each element gets saved to the work space as you adjust it. If you want to save the elements current configuration (size, position, text, name) as a template, the option to save as template is at the bottom of the properties. You can also find the delete element option there.

2.4 Preview and Save

Once you have completed your elements and survey settings, you can preview your survey in the Preview and Save bar. There is also an option to clear the survey, rename the survey, or load from an existing survey template. There are then two save options;

Save to File- will store a file version of this survey slide that can be loaded again later for editing but WILL NOT save the survey slide to your stimulus presentation flow of the study you are working in.

Save to Study- will add this current survey slide design to the current study's stimulus presentation.

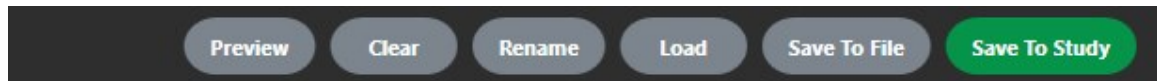


Note that the current survey slide name is shown in the bottom right corner of the interface.

3. Element Property Details


3.1 Rich Text Message

The rich text message editor is a fully flexible text editor so you can control the font, size, color, and content of the message you would like to display to the respondent. Because the Rich Text Message opens the 'Edit Rich Text' box, there are no additional properties in the properties section for this element.

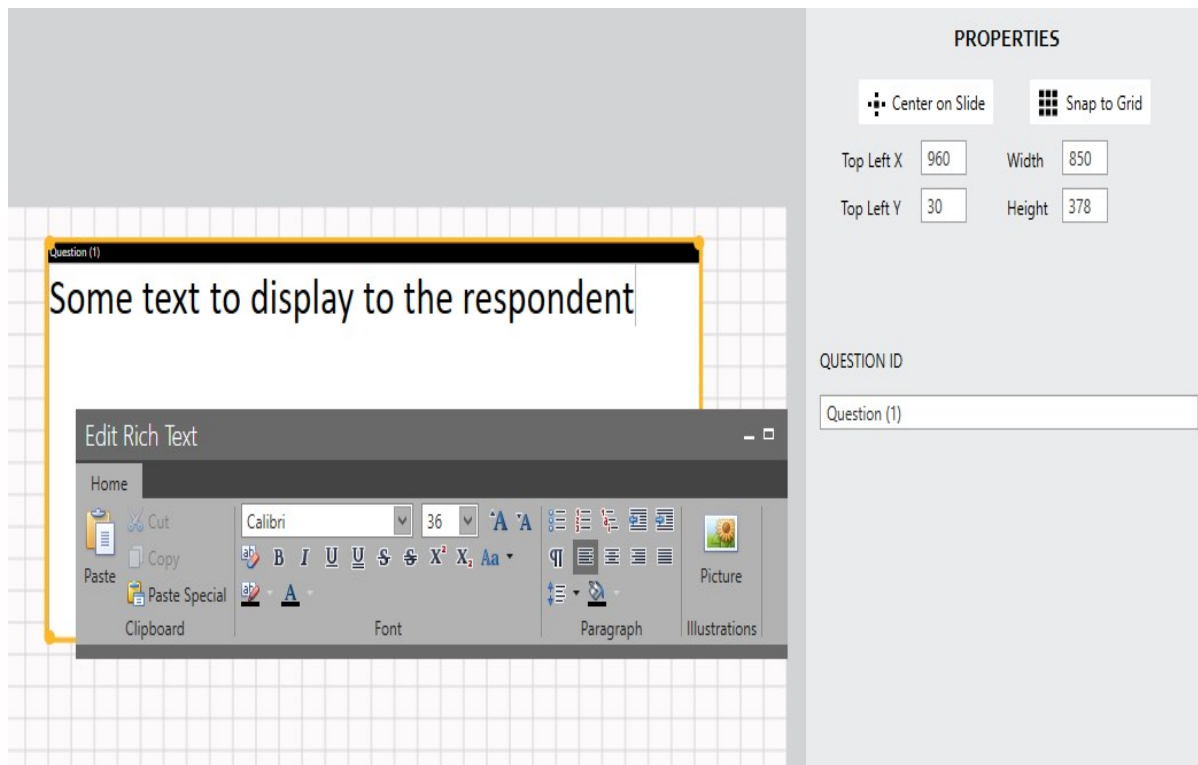


3.2 Plain Text Message

The Plain Text Message box is very similar to the Rich Text Message box but without the 'Edit Rich Text' box. This could preferable over the rich text message to quickly add a message element but also serves the purpose of displaying a message to the respondent. The text is written in the 'Message' box in the Properties section where the text size can also be adjusted.

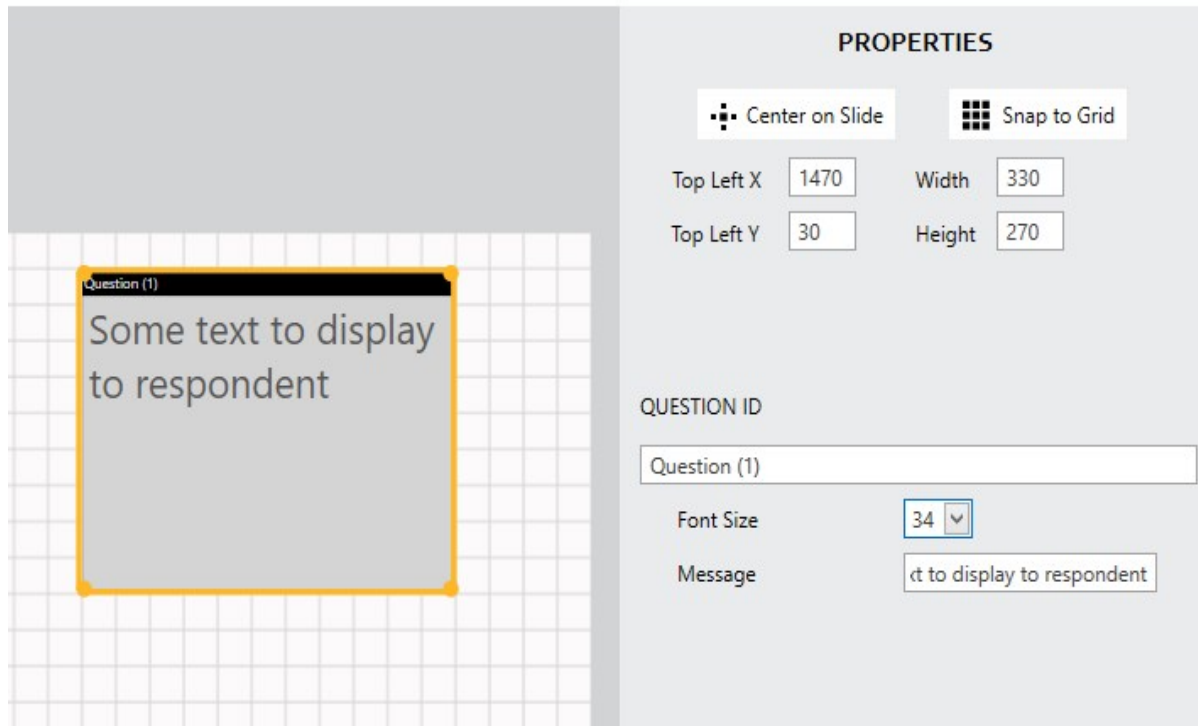

3.3 Image Element

Selecting the Image element will open the file explorer and allow you choose a .jpeg image file. You can add a header and adjust it's font size in the properties section. Note that you can also swap the image file location in the properties where it says Image File.

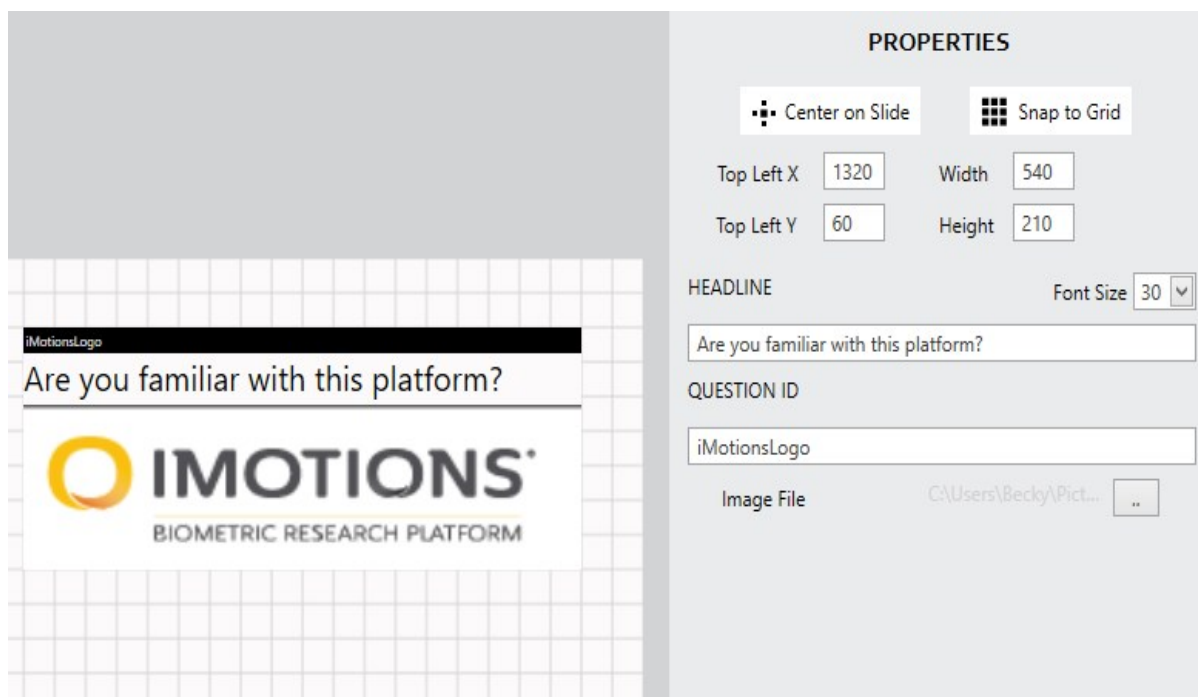

3.4 Scale

The sliding scale element allows the respondent to slide to either discrete values or on a continuous bar. You can adjust the labels that are shown to the respondent to represent the regions of the bar as well as the value each selection is worth.

The image shows a survey design interface with a grid background. On the left, a preview of a 'Scale Question' is shown. The question text is 'Slide to select'. Below it is a horizontal scale bar with five discrete points labeled -2, -1, 0, 1, and 2. Under these labels are the corresponding response options: 'Definitely No' for -2, 'Somewhat No' for -1, 'I'm Not Sure' for 0, 'Somewhat Yes' for 1, and 'Definitely Yes' for 2. A small slider icon is positioned at the value 1. On the right, a 'PROPERTIES' panel is visible. It contains settings for the question, including 'Center on Slide' and 'Snap to Grid' checkboxes, 'Top Left X' (690), 'Width' (1200), 'Top Left Y' (150), and 'Height' (240). The 'HEADLINE' is 'Slide to select' with a 'Font Size' of 30. The 'QUESTION ID' is 'Scale Question' with a 'Font Size' of 20. Other settings include 'Required' (No), 'Scale Type' (Discrete), 'Show Values' (Yes), 'Show Labels' (Yes), 'Default Value' (Somewhat Yes), and 'Number of bins' (10). There are 'Edit' buttons for 'Labels & values' and 'Histogram bins'.


The exported survey results you generate after running your study will display the selected value and label the respondent chose for this Question ID so make sure you design these accordingly. Values can be used to weight certain responses. For example, here the two negative answers were made to be negative numbers.

If you choose discrete, it will automatically snap the respondents selection to one of the values. Choosing "continuous" allows you to setup histogram bins so that you can define the granularity you are interested in viewing the respondents' choices. For example, here there are 10 bins. The respondents selection will fit within one of these bins so you can have a distribution of responses across 10 bins in the exported results.


3.5 Single / Multiple Choice

A common method of surveying is to ask the respondent to select between choices. You can set up one or more multiple choice questions in this one element so that they are all aligned perfectly.

Open the Questions editor to choose how many and what questions to ask the respondent.

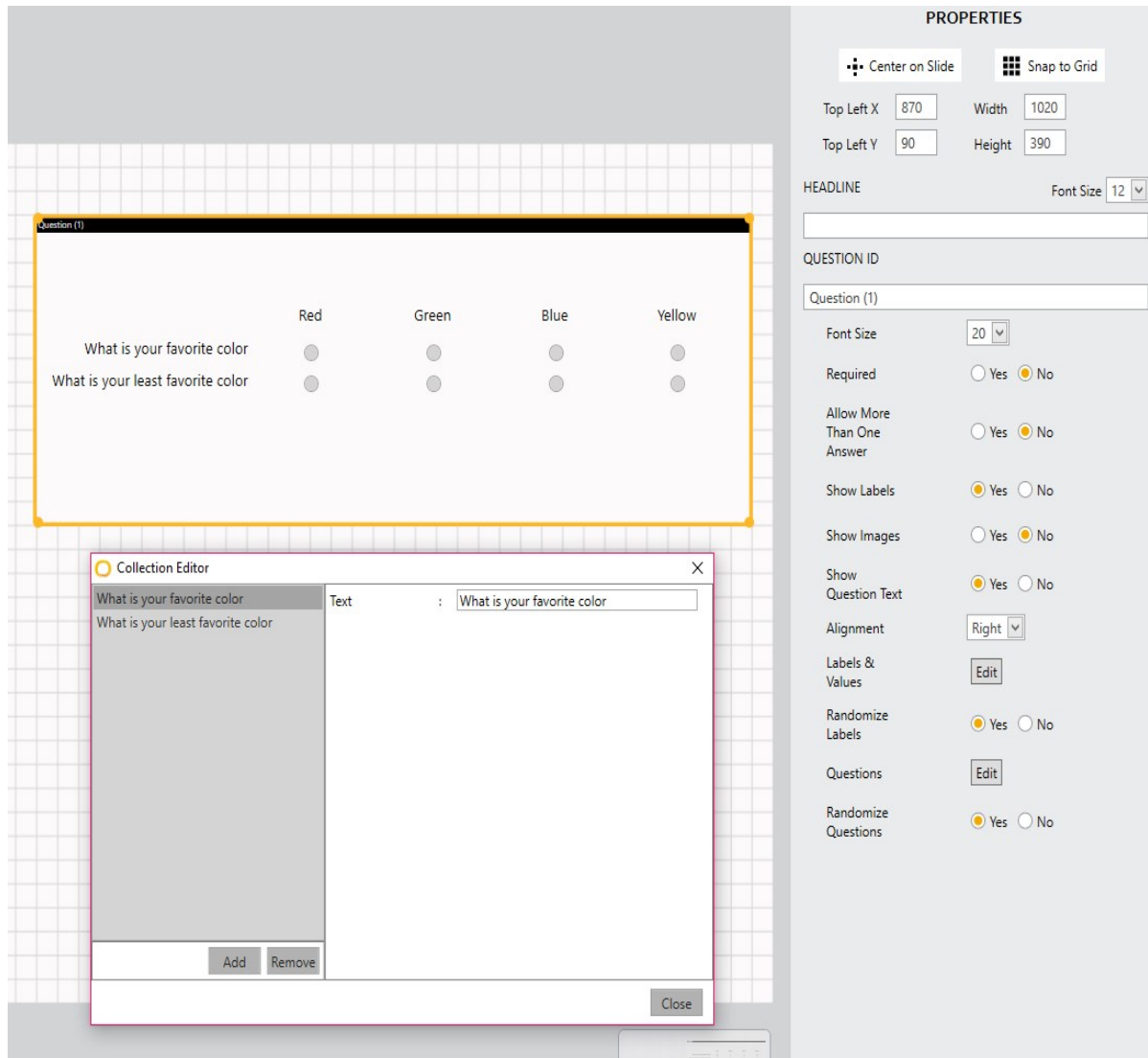

Open the labels and values editor to adjust what response options the respondent will see (labels) and what each selection is worth (values).

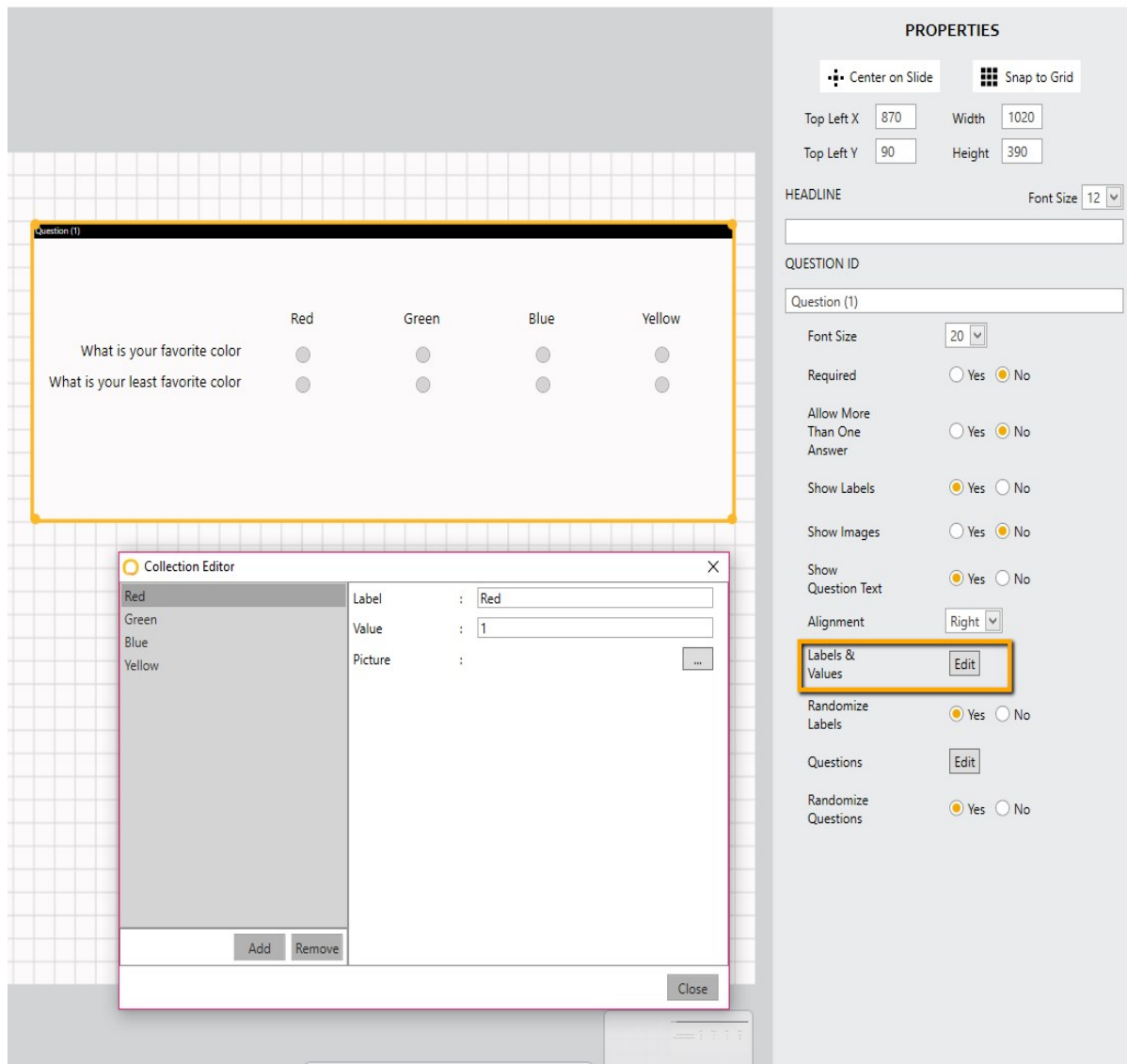

Other properties include: the ability to require a response to all questions in this element before advancing (Required), the ability to choose multiple options to each question (Allow More than One Answer), show the labels of response options or not (Show Labels), show the questions or not (Show Question Text), align the question text (Alignment), and randomize the questions in the element.


Required	☐ Yes ☒ No
Allow More Than One Answer	☐ Yes ☒ No
Show Labels	☒ Yes ☐ No
Show Images	☐ Yes ☒ No
Show Question Text	☒ Yes ☐ No
Alignment	Right v
Labels & Values	Edit
Randomize Labels	☒ Yes ☐ No
Questions	Edit
Randomize Questions	☒ Yes ☐ No

NOTE ABOUT RANDOMIZATION AND REPLAYS / EYE TRACKING

If you choose to randomize either the questions or the labels, while the respondents will view the randomization appropriately and the responses are stored in association with the correct question, all video replays in iMotions will be associated with one image (one question orientation). Keep this in mind if you are drawing any areas of interest or performing eye tracking analysis. Please discuss with a product specialist if you have concerns about this.

3.6 Vertical Selection

The vertical selection element is very similar to the single / multiple choice element but displays the questions across the top and the labels down the side. Please see the above section for a description of the properties.

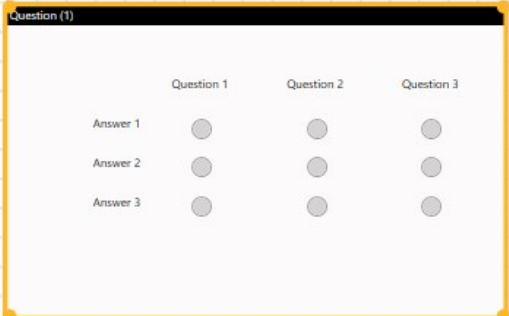

PROPERTIES

☒ Center on Slide

☒ Snap to Grid

Top Left X

Width

Top Left Y

Height

HEADLINE Font Size

QUESTION ID

Font Size

Required ☐ Yes ☒ No

Allow More Than One Answer ☐ Yes ☒ No

Show Labels ☒ Yes ☐ No

Show Images ☒ Yes ☐ No

Show Question Text ☒ Yes ☐ No

Alignment

Labels & Values


Randomize Labels ☐ Yes ☒ No

Questions

Randomize Questions ☐ Yes ☒ No

3.7 Text Response

You can ask the respondent to write a short answer free response in a text response element. The Headline property is one way you could write a prompt for the respondent. Alternatively, leave this blank and add a Rich Text Element to the slide as well. Additional properties include the option to require a response and the line length which determines how far their text will extend before defaulting to the next line.

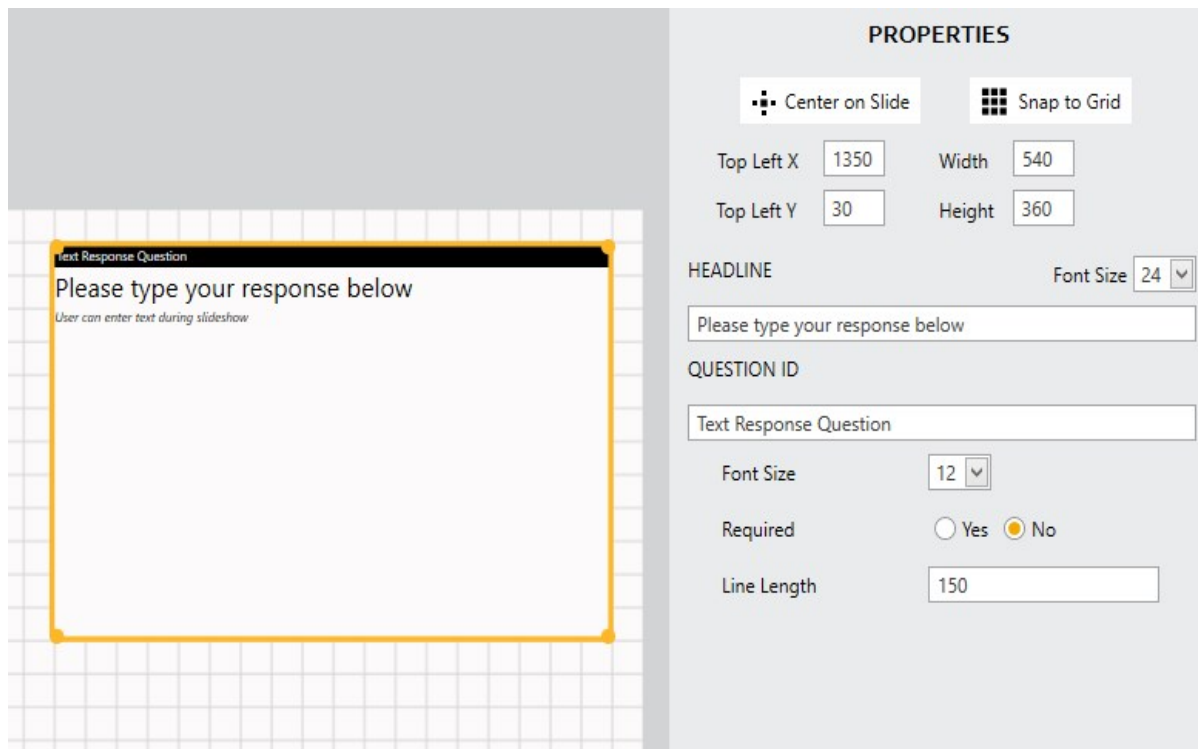


3.8 Key Button Response

The Key Button Response element allows you to setup key presses the respondent will use to respond to a question with. In this case, the labels and values editor will ask for a keyboard option the respondent can answer with in addition to the question label as it will appear in the exported data and it's associated value.


If you select Advance Upon Keypress, iMotions will automatically progress to the next stimulus when it receives a key press that is in the configured list for this element.

Key Button Response questions cannot be combined on the same slide with other forms of response elements.

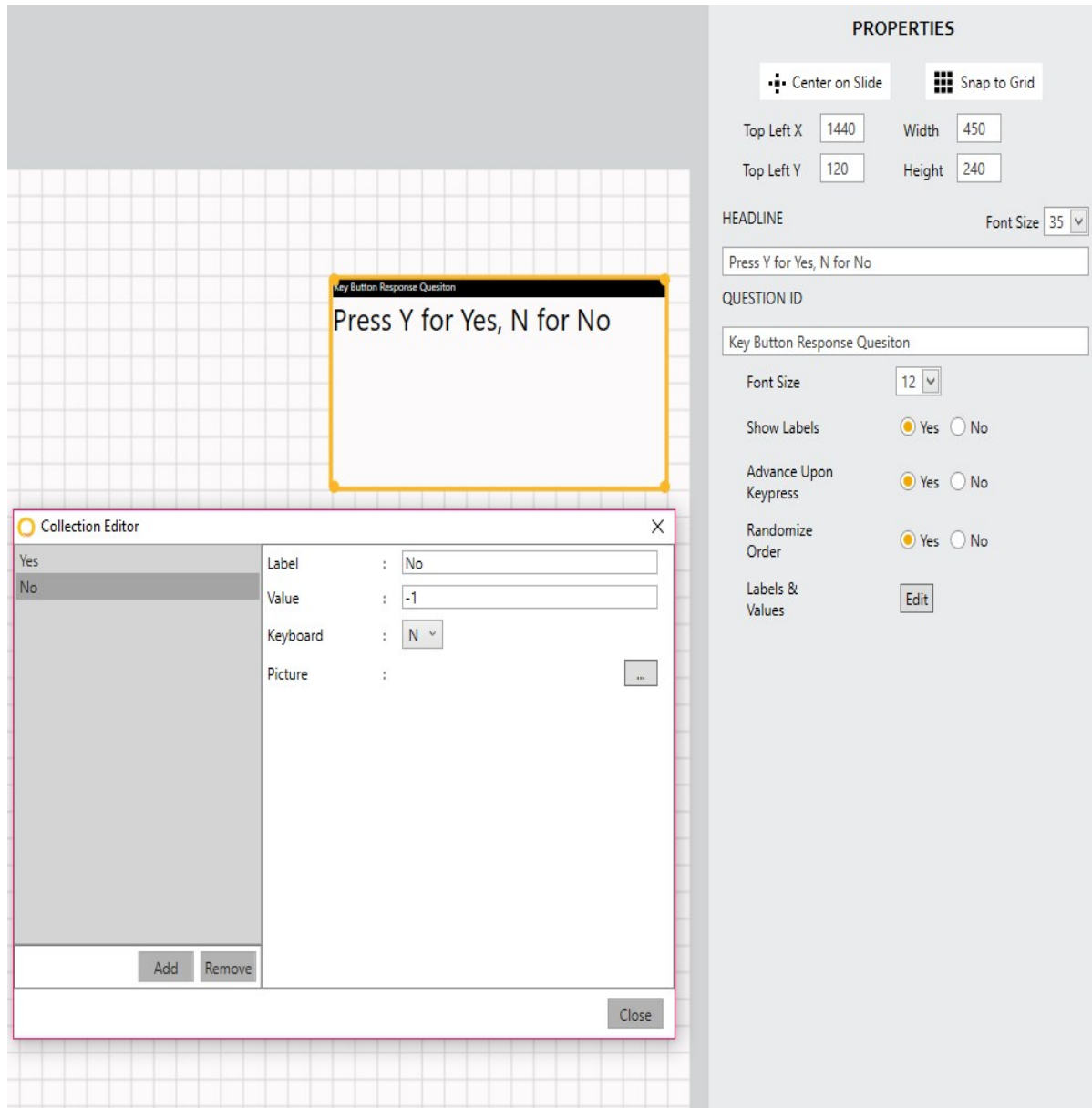

3.9 Survey Layer controller


In iMotions 7.1 its possible to utilize the survey layer controller, meaning you have the option to bring a object to the front position in your survey. Pre 7.1, the latest created object would occupy the front position.

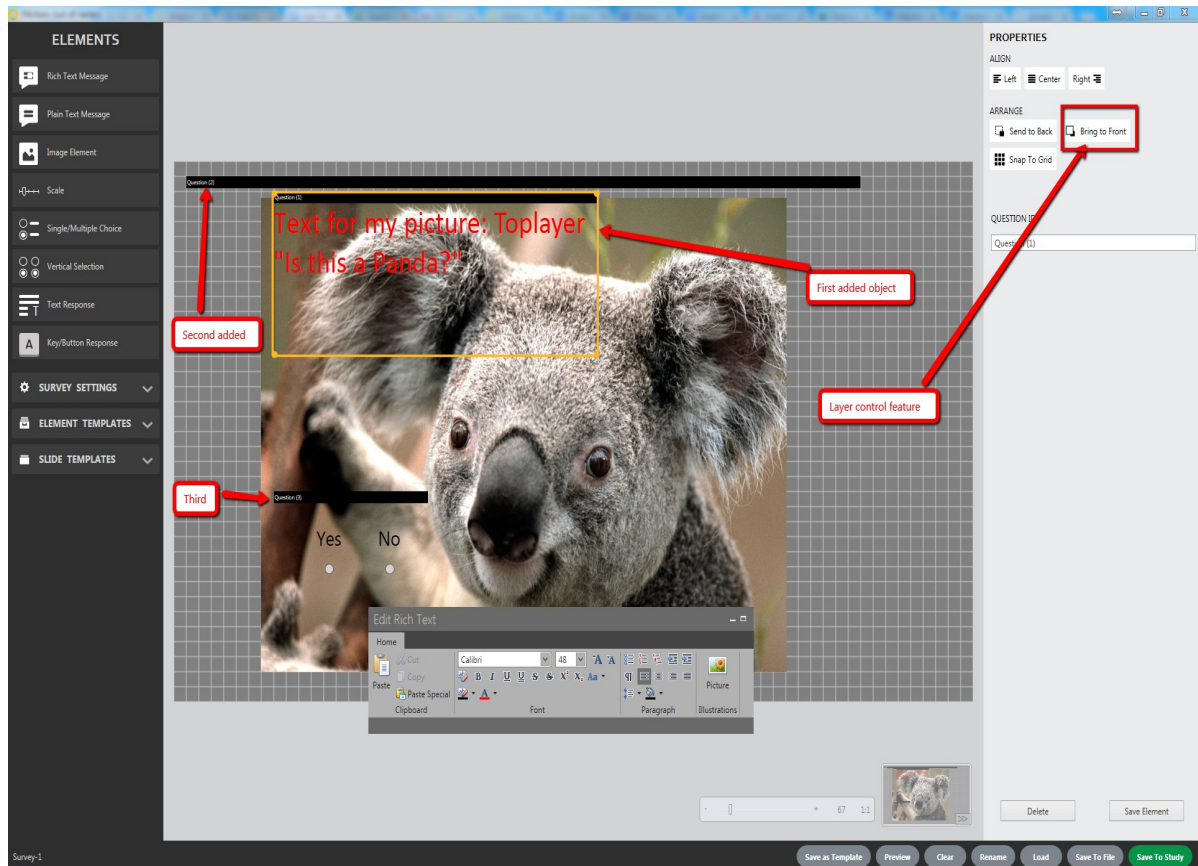

In the example above, the Text field was added first, then the picture, then the Question field. However the button "Bring to front" was used on first, the Text field, then the Question field. So currently the Question field is "Layer 1" but since its not overlapping the Text field, this doesnt affect the layout.


4. Survey Settings

The survey settings allow you to adjust the background color of the slide, the font color, resolution, whether to display the grids on the work space while you work, whether to auto align the elements to the grid as you move them, and to adjust the grid size.

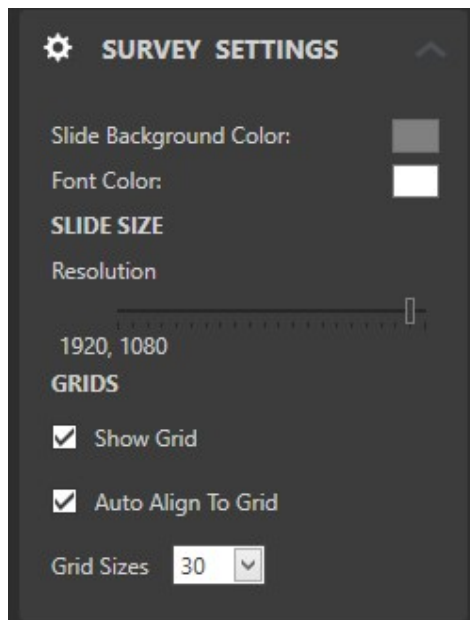

5. Templates

There are two groups of templates: Elements and Slides. The Element templates will start out blank, but you can populate it by creating an element and setting up it's properties then choosing 'save as template' in the properties section for the selected element. You will then choose a location to save the template and if you keep the default location, will see the element here.

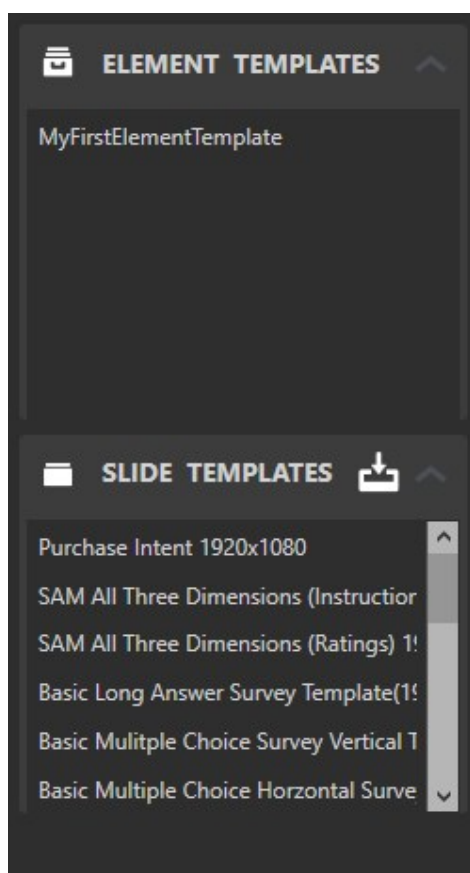

The slide templates comes populated with some default survey slides. You can learn more about these in the survey section of the help center. Note that selecting a template will clear the contents of the work space. If you want to add additional survey slides to this section, choose the 'save to file' option in the preview and save section of the interface.

You can learn more about the included templates here:

[Survey Section of the Help Center](#)

