

Table of Contents

iMotions Lab > Data Collection > Importing Data into iMotions > Eye Tracking Glasses

Importing Eye Tracking Glasses Recordings/Files	2
---	---

Importing Eye Tracking Glasses Recordings/Files

iMotions offers a range of import options for eye-tracking glasses manufactured by iMotions' hardware partners. Suppose you wish to use a device that is not officially supported by iMotions. In that case, you have the option to convert your recorded data into a format that is recognized by iMotions. This will then allow you to import your data into iMotions and use all of iMotions analysis capabilities nonetheless.

This Help Center article gives information about the

- [Expected file format](#)
- [Expected folder structure](#)
- [How to post import eye tracking glasses data into iMotions](#)
- [Next steps](#)

Expected file format

For every respondent, you need

- one video recording from the glasses' scene camera with the extensions `.mp4`, `.avi`, or `.wmv`
- one `.csv` file with the following eye-tracking data
 - *Timestamp* (continuous timestamp in milliseconds)
 - *GazeLeftX* (left horizontal gaze coordinates in pixels. Origin must be the video's upper left corner, value range must equal the resolution of the video)
 - *GazeLeftY* (left vertical gaze coordinate in pixels. Origin must be the video's upper left corner, value range must equal the resolution of the video)
 - *PupilDiaLeft* (diameter of the left pupil in millimeters)
 - *GazeRightX* (right horizontal gaze coordinates in pixels. Origin must be the video's upper left corner, value range must equal the resolution of the video)
 - *GazeRightY* (right vertical gaze coordinates in pixels. Origin must be the video's upper left corner, value range must equal the resolution of the video)

- *PupilDiaRight* (diameter of the right pupil in millimeters)

All of the fields are required and the naming convention as well as the units are obligatory. You can use the attached template.

In case your device does not differentiate between data from the left and right eye we recommend duplicating the data.


In case your device does not provide pupil diameter data, we recommend filling in an arbitrary constant value instead.

iMotions will not typically provide a converter for your data; it is your responsibility to correctly format your recordings to the format that iMotions can read.

Expected folder structure

For every respondent, save the .csv file as well as the video file in one folder per respondent. The name of this folder will later be used as the Respondent Name in iMotions.


All these respondent folders need to then be saved into one study folder.

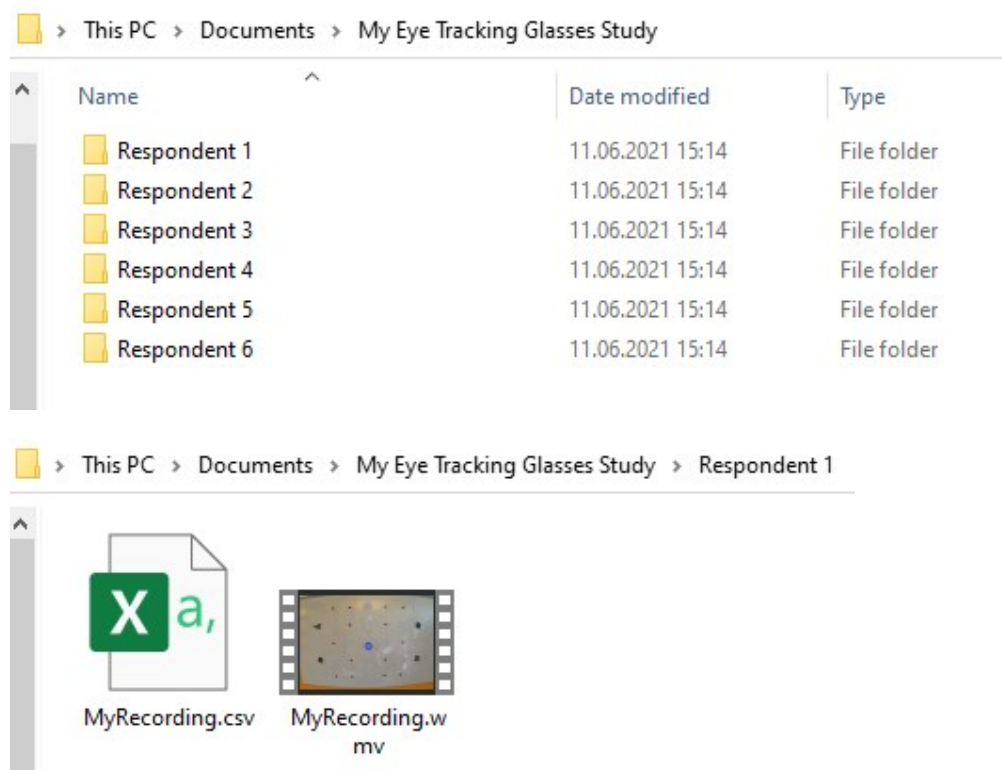

How to post import eye tracking glasses data into iMotions

For iMotions 10 (or newer), Eye Tracking Glasses data are imported through the New Study interface.

1. Create a New Study in iMotions
2. Select Glasses study and click Next
3. Select Import glasses data and click Next
4. Select your Eye tracking glasses/model
5. Choose to import files as a New study or to Merge into an existing compatible study (with the same glasses, settings, etc.)
6. Select the Data folder that contains the the respondent folders and files (described above). Confirm the folder structure complies with the description for the selected eye tracking glasses.
7. Click Import to create respondents with eye tracking glasses data/video recordings in iMotions for each respondent folder

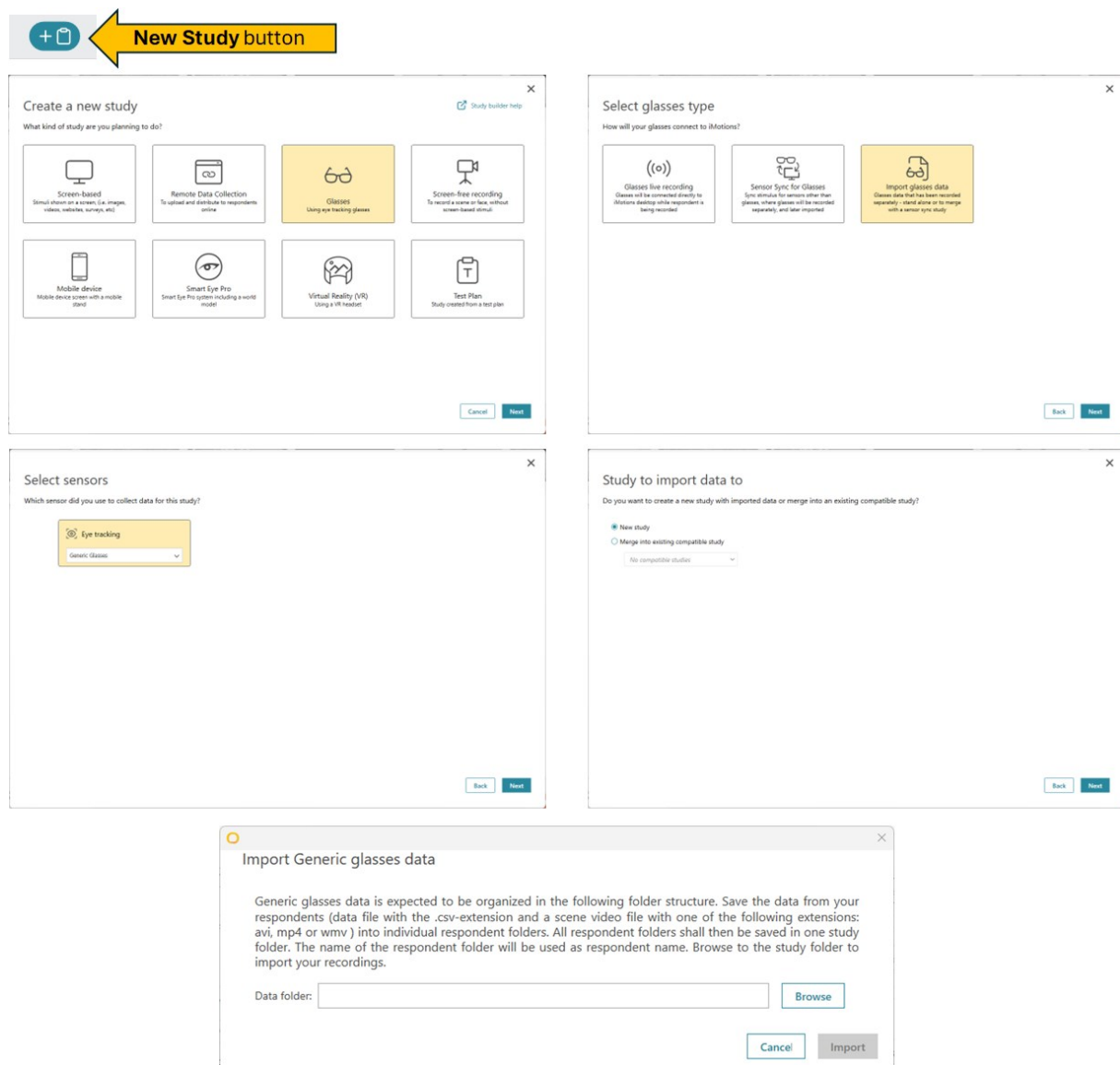

Next steps

Once the data is imported, it is available for all eye-tracking analytics such as:

- Gaze Mapping
- Areas of Interest
- Fixation classification

and synchronization with other sensors.

